

# Macroresults through Microarrays

John C. Rockett, Reproductive Toxicology Division (MD-72), National Health and Environmental Effects Research Laboratory, Office of Research and Development, US Environmental Protection Agency, Research Triangle Park, 2525 East Highway 54, Durham, NC 27711, USA;  
tel: +1 919 541 2071, fax: +1 919 541 4017, e-mail: rockett.john@epa.gov

The third enactment of Cambridge Healthtech Institute's *Macroresults through Microarrays* meeting was held in Boston (MA, USA) from 29 April–1 May 2002. The subtheme of this year's meeting was 'advancing drug discovery', a widely touted application for array technology.

## The evolution of microarrays

If you were asked 'Who first conceived of the idea of microarrays', who would come to mind? Mark Schena perhaps, first author of the seminal 1995 paper on cDNA arrays [1]? Maybe Pat Brown, Schena's then supervisor? Or perhaps Stephen Fodor, the primary driver behind Affymetrix's (<http://www.affymetrix.com>) oligonucleotide-based platform [2]. Brits might even chant the name of Ed Southern [3]. Well, according to Roger Ekins (University College London Medical School; <http://www.ucl.ac.uk/medicine/>) all these answers would be wrong. It was in fact Ekins and his colleagues who first conceived of and patented 'a new generation of ultrasensitive, miniaturized assays for protein and DNA–RNA measurement based on the use of microarrays' in the mid 1980s [4]. The concept and potential of array technology was more fully described in a later publication, in which Ekins *et al.* [5] concluded that antibody microspots of  $\sim 50 \mu\text{m}^2$  could be achieved, and that as many as 2 million different immunoassays could, in principle, be accommodated on a surface area of  $1 \text{ cm}^2$ .

## Technological innovation

In practice, it took a different biological molecule (DNA), a different research

group, and a leap into microfabrication technology to even begin approaching these kinds of densities [Affymetrix patent 6045996 talks of one million spots  $\text{cm}^{-2}$ ]. Of course, advancing technology is one of the driving engines behind the genomics juggernaut, and we are already seeing '4th generation' machines for fabricating DNA chips. If the company representatives at this meeting are to be believed (and their cases seemed strong), spotting is out, and *in situ* fabrication of oligonucleotide-based 'iterative custom arrays' is in. Whether you go with the Combimatrix's (<http://www.combimatrix.com>) electrochemically directed synthesis and detection system, febit's (<http://www.febit.com>) Geniom® technology, or Nimblegen's (<http://www.nimblegen.com>) Maskless Array Synthesizer technology is a matter of personal choice. However, each of these machines provides the flexibility to design variable length oligonucleotide probes from sequences inputted by the user, and then perform *in situ* synthesis of an array. Each system also boasts unique advantages. For example, Combimatrix's biological array processor is a semiconductor coated with a 3D layer of porous material in which DNA, RNA, peptides or small molecules can be synthesized or immobilized within discrete test sites, while febit's Geniom One® is a fully integrated gene-expression analysis system with minimal user hands-on time – the probe sequences are programmed, the RNA samples inserted, and the gene expression data is pumped out a few hours later.

## Cell- and tissue-based arrays

Array technology is in most people's minds firmly linked with gene-expression profiling. Fewer are aware that cell- and tissue-based arrays have been developed, and how they can provide a vital extra dimension to research. In support of this, Barry Bochner gave an update on the cell-based array system that Biolog (<http://www.biolog.com>) has produced for simultaneously measuring the effects of one gene in the cell under thousands of growth conditions (see [6] for further details). David Walt (Tufts University; <http://www.tufts.edu/>) is developing single live cell arrays using optical imaging fiber (OIF) technology. An array of microwells is fabricated on the face of an OIF at densities of up to 10 million wells  $\text{cm}^{-2}$ . Cells are then added to the wells and disperse at an average of one cell per well. Physiological and genetic responses of each cell are measured via fluorescence produced by reporter genes (e.g. *lacZ*, *gfp*). Assays performed so far include yeast live or dead cell assay, microenvironment pH and  $\text{O}_2$  measurements, promotor responses using the *lacZ* and *phoA* reporter genes, and protein–protein interactions using the yeast two-hybrid system. The main advantage of this system is that the cells remain alive during the assay, which means a real-time timecourse can be performed and/or the array passed from sample to sample. This would be useful in, for example, the scanning of a combinatorial drug library for specific physiological effects.

Tissue arrays are a useful complementary technology to DNA arrays because they can be used to help validate and

understand the biological and medical significance of gene changes discovered using standard DNA arrays. For example, an array of tumor tissues can be screened for the protein (using immunohistochemistry), message (using *in situ* hybridization) and copy number (using comparative genomic hybridization) of a gene of interest, to determine if expression of the gene (or lack thereof) is related in any way to survival. They can also be used to predict the probability of clinical failure of lead compounds as a result of toxicity by evaluating the distribution of the drug targets in normal tissue. Spyro Mousses and his co-workers at the National Human Genome Research Institute (<http://www.nhgri.nih.gov/index.html>) have built such arrays, including a multi-tumor array (~5000 specimens, and sections from 36 normal and 800 metastatic tissues) and a normal tissue array (76 tissue and 332 cell types).

### The problem with proteins

It has been said that genomics tells us what might happen, transcriptomics indicates what should happen, and proteomics shows what is happening. The impact of functional proteomics on pharmaceutical R&D is rapidly increasing, and protein arrays are being used increasingly in both basic and applied research. Their use lies not only in comparative protein expression and interaction profiling, but also in diagnostics and drug discovery. However, an increasing number of researchers have found that protein arrays, like their cousins the DNA arrays, present several practical obstacles relating to their production and use. For example, in using *Escherichia coli* to produce recombinant eukaryotic proteins from a single expression vector, multiple protein products are often produced, suggesting mixes of truncated or otherwise altered proteins. There is also the obvious concern that the proteins might not be modified in a similar manner to

eukaryotic systems. Also, an optimal method for depositing and binding proteins to the selected substrate is yet to be determined, as is the best way to ensure that they are bound in a correctly folded, active conformation.

Several companies have been addressing these problems. Prolinx (<http://www.prolinxinc.com>) is one such company, and Karin Hughes described their Versalinx™ chemistry for producing protein, peptide and small-molecule arrays. Versalinx™ uses solution-phase conjugation followed by immobilization, resulting in functional orientation of proteins and peptides on the substrate surface. It also offers the valuable additional benefit of exhibiting low non-specific binding. Sense Proteomic (<http://www.senseproteomic.com>) is also among those addressing these problems to develop robust protein arrays for drug discovery and clinical applications and has developed functional protein array formats based on specific disease tissues. Subtractive hybridization is used to identify genes with altered expression in breast tumor and cystic fibrosis compared to normal tissue. A high throughput cloning strategy (COVET™) is then used to produce libraries of genes that are tagged, cloned, expressed, purified and finally immobilized on glass slides. Initial validation studies have shown that the vast majority of the immobilized proteins do indeed retain biological function.

Stefan Schmidt and his company (GPC Biotech; <http://www.gpcbiotech.de>) have moved past the platform development stage and, with their focus firmly on drug discovery, are currently developing kinase-profiling arrays. Kinases are important targets for pharmaceutical drug discovery and therapy, and GPC's aim is to simultaneously detect multiple kinases, obtain activity profiles for different cell types, or analyze the ability of drug candidates to inhibit kinase activity. To do this, recombinant kinase substrates are immobilized on

membranes, incubated with purified kinase, and the substrates measured for the degree of phosphorylation.

### Summary

Meetings like this, packed with exciting discoveries and intriguing and interesting innovation, heavily emphasize the pace at which biotechnology is advancing, to the extent that the number of options for genomic and proteomic researchers can become overwhelming. Although data analysis is perhaps the greatest current concern for array users, an increasing challenge will be to determine the approaches and technology that really work, and to do it in a timely manner.

### References

- 1 Schena, M. *et al.* (1995) Quantitative monitoring of gene expression patterns with a complementary DNA microarray. *Science* 270, 467-470
- 2 Fodor, S.P. *et al.* (1991) Light-directed, spatially addressable parallel chemical synthesis. *Science* 251, 767-773
- 3 Southern, E.M. *et al.* (1992) Analyzing and comparing nucleic acid sequences by hybridization to arrays of oligonucleotides: evaluation using experimental models. *Genomics* 13, 1008-1017
- 4 Ekins, R.P. (1987) US Patent Application 8 803 000
- 5 Ekins, R. *et al.* (1989) High specific activity chemiluminescent and fluorescent markers: their potential application to high sensitivity and 'multi-analyte' immunoassays. *J. Biolum. Chemilum.* 4, 59-78
- 6 Rockett, J.C. (2002) Chip, chip, array! Three chips for post-genomic research. *Drug Discov. Today* 7, 458-459

### Acknowledgements

I would like to thank Mary Ann Brown (Cambridge Healthtech Institute) and David Dix (US EPA) for critical review of this manuscript prior to submission. This document has been reviewed in accordance with US Environmental Protection Agency policy and approved for publication. Mention of companies, trade names or products does not signify endorsement of such by the EPA.

N. Leigh Anderson  
Ricardo Esquer-Blasco  
Jean-Paul Hofmann  
Norman G. Anderson

Large Scale Biology Corporation,  
Rockville, MD

## A two-dimensional gel database of rat liver proteins useful in gene regulation and drug effects studies

A standard two-dimensional (2-D) protein map of Fischer 344 rat liver (F344MST3) is presented, with a tabular listing of more than 1200 protein species. Sodium dodecyl sulfate (SDS) molecular mass and isoelectric point have been established, based on positions of numerous internal standards. This map has been used to connect and compare hundreds of 2-D gels of rat liver samples from a variety of studies, and forms the nucleus of an expanding database describing rat liver proteins and their regulation by various drugs and toxic agents. An example of such a study, involving regulation of cholesterol synthesis by cholesterol-lowering drugs and a high-cholesterol diet, is presented. Since the map has been obtained with a widely used and highly reproducible 2-D gel system (the Iso-Dalt® system), it can be directly related to an expanding body of work in other laboratories.

Lal et al., 09/002,485, filed December 31, 1997  
(PF-0459)

Exhibit "L" attached to Declaration of John C.  
Rockett, Ph.D.

### Contents

1 Introduction.....	907
2 Material and methods.....	908
2.1 Sample preparation.....	908
2.2 Two-dimensional electrophoresis.....	909
2.3 Staining.....	909
2.4 Positional standardization.....	909
2.5 Computer analysis.....	909
2.6 Graphical data output.....	910
2.7 Experiment LSBC04.....	910
3 Results and discussion.....	910
3.1 The rat liver protein 2-D map.....	910
3.2 Carbamylated charge standards computed pI's and molecular mass standardization.....	911
3.3 An example of rat liver gene regulation: Chol- esterol metabolism.....	911
3.3.1 MSN 413 (putative cytosolic HMG-CoA synthase) and sets of spots regulated co- ordinately or inversely.....	911
3.3.2 MSN 235 and coregulated spots.....	912
3.3.3 An example of an anti-synergistic effect.....	912
3.3.4 Complexity of the cholesterol synthesis pathway.....	912
4 Conclusions.....	912
5 References.....	912
6 Addendum 1: Figures 1-13.....	914
7 Addendum 2: Tables 1-4.....	923
Table 1. Master table of proteins in rat liver data- base.....	923
Table 2. Table of some identified proteins.....	928
Table 3. Computed pI's of two sets of carbamylated protein standards: rabbit muscle CPK and human Hb.....	929
Table 4. Computed pI's of some known proteins re- lated to measured CPK pI's.....	930

### 1 Introduction

High-resolution two-dimensional electrophoresis of proteins, introduced in 1975 by O'Farrell and others [1-4], has been used over the ensuing 16 years to examine a wide variety of biological systems, the results appearing in more than 5000 published papers. With the advent of computerized systems for analyzing two-dimensional (2-D) gel images and constructing spot databases, it is also possible to plan and assemble integrated bodies of information describing the appearance and regulation of thousands of protein gene products [5, 6]. Creating such databases involves amassing and organizing quantitative data from thousands of 2-D gels, and requires a substantial commitment in technology and resources.

Given the long-term effort required to develop a protein database, the choice of a biological system takes on considerable importance. While *in vitro* systems are ideal for answering many experimental questions, especially in cancer research and genetics, our experience with cell cultures and tissue samples suggests that some *in vivo* approaches could have major advantages. In particular, we have noticed that liver tissue samples from rats and mice appear to show greater quantitative reproducibility (in terms of individual protein expression) than replicate cell cultures. This is perhaps a natural result of the homeostasis maintained in a complete animal vs. the well-known variability of cell cultures, the latter due principally to differences in reagents (e.g., fetal bovine serum), conditions (e.g., pH) and genetic "evolution" of cell lines while in culture. It is also more difficult to generate adequate amounts of protein from cell culture systems (particularly with attached cells), forcing the investigator to resort to radioisotope-based or silver-based stain-detection methods. While these methods are more sensitive (sometimes much more sensitive) than the Coomassie Brilliant Blue (CBB) stain typically used for protein detection in "large" protein samples, they are generally more variable, more labor-intensive and, in the case of radiographic methods, may generate highly "noisy" images, due to the properties of the films used. By contrast, large protein samples can easily be prepared from liver using urea/Nonidet P-40 (NP-40) solubilization and stained with CBB, which has the advantage of being easily reproducible [8]. Finally, there remains the question of the "truthfulness" of many *in vitro* systems as compared to their *in vivo* analogs; how great are the changes caused by the introduction into a cul-

Correspondence: Dr. N. Leigh Anderson, Large Scale Biology Corporation, 9620 Medical Center Drive, Rockville, MD 20850, USA

Abbreviations: CBB, Coomassie Brilliant Blue; CPK, creatine phosphokinase; 2-D, two-dimensional; IEF, isoelectric focusing; MSN, master spot number; NP-40, Nonidet P-40; SDS, sodium dodecyl sulfate

ture and the associated shift to strong selection for growth, and how do these affect experimental outcomes? Hence the apparent advantages of *in vitro* systems, in terms of experimental manipulation, may be counterbalanced by other factors relating to 2-D data quality.

There is a second important class of reasons for exploring the use of an *in vivo* biological system such as the liver. Historically, there have been two broad approaches to the mechanistic dissection of biochemical processes in intact cellular systems: genetics (a search for informative mutants) and the use of chemical agents (drugs and chemical toxins). Both approaches help us to understand complex systems by disrupting some specific functional element and showing us the result. With the development of techniques for genetic manipulation and cloning, the genetic approach can be effectively applied either *in vitro* or *in vivo*, although the *in vitro* route is usually quicker. The chemical approach can also be applied to either sort of biological system; here, however, the bulk of consistently acquired information is in experimental animals (rats and mice). While most biologists know a short list of compounds having specific, experimentally useful effects (e.g., inhibitors of protein synthesis, ionophores, polymerase inhibitors, channel blockers, nucleotide analogs, and compounds affecting polymerization of cytoskeletal proteins), there is a much larger number of interesting chemically-induced effects, most of them characterized by toxicologists and pharmacologists in rodent systems. Just as a thorough genetic analysis would involve saturating a genome with mutations, it is possible to imagine a saturating number of drugs, the analysis of whose actions would reveal the complete biochemistry of the cell. While organized drug discovery efforts usually target specific desired effects, the nature of the process, with its dependence on screening large numbers of compounds, necessarily produces many unanticipated effects. It is therefore reasonable to suppose that the required broad range of compounds necessary to achieve "biochemical saturation" may be forthcoming; in fact, it may already exist among the hundreds of thousands of compounds that failed to qualify as drugs.

Among organs, the liver is an obvious choice for the study of chemical effects because of its well-known plasticity and responsiveness. The brain appears to be quite plastic (e.g. [7]), but it is a complicated mixture of cell types requiring skillful dissection for most experiments. The kidney, while quite responsive, also presents a potentially confounding mixture of cell types. The liver, by contrast, is made up of one predominant cell type which is easy to solubilize: the hepatocyte, representing more than 95% of its mass. Most importantly, the liver performs many homeostatic functions that require rapid modulation of gene expression. It appears that most chemical agents tested affect gene expression in the liver at some dosage (N. Leigh Anderson, unpublished observations), an interesting contrast to our earlier work with lymphocytes, for example, which seem to be much less responsive. Such results conform to the expectation that cells with a homeostatic, physiological role should be more plastic than cells differentiated for a purpose dependent on the action of a limited number of specific genes.

The liver also allows the parallels between *in vitro* and *in vivo* systems to be examined in detail. Significant progress

has been made in the development of mouse, rat and human hepatocyte culture systems, as well as in precision-cut tissue slices. Using such an array of techniques, it is possible to assemble a matrix of mammalian systems including mouse and rat *in vivo* on one level and mouse, rat and human *in vitro* on a second level, and to compare effects between species and between systems. This approach allows us to draw informed conclusions regarding the biochemical "universality" of biological responses among the mammals, and to offer some insight into the validity of *in vitro* approaches for toxicological screening. We believe this data will be necessary if *in vitro* alternatives are to achieve wide usage in government-mandated safety testing of drugs, consumer products and industrial and agricultural chemicals.

A number of interesting studies have been published using 2-D mapping to examine effects in the rodent liver. A number of investigators have made use of the technique to screen for existing genetic variants [8-11] or induced mutations [12-14], mainly in the mouse. This work builds on the wealth of genetic information available on the mouse and its established position as a mammalian mutation-detection system. While some studies of chemical effects have been undertaken in the mouse [15-17], most have used the rat [18-23]. The examination of the cytochrome p-450 system, in particular, has been carried out almost exclusively on the rat [24, 25].

These considerations lead us to conclude that rodent liver offers the best opportunity to systematically examine an array of gene regulation systems, and ultimately to build a predictive model of large-scale mammalian gene control. The basic underlying foundation of such a project is a reliable, reproducible master 2-D pattern of liver, to which ongoing experimental results can be referred. In this paper, we report such a master pattern for the acidic and neutral proteins of rat liver (pattern F344MST3). In future, this master will be supplemented by maps of basic proteins, and analogous maps of mouse and human liver.

## 2 Materials and methods

### 2.1 Sample preparation

Liver is an ideal sample material for most biochemical studies, including 2-D analysis. A sample is taken of approximately 0.5 g of tissue from the apical end of the left lobe of the liver. Solubilization is effected as rapidly as practical; a delay of 5-15 min appears to cause no major alteration in liver protein composition if the liver pieces are kept cold (e.g., on ice) in the interim. In the solubilization process, the liver sample is weighed, placed in a glass homogenizer (e.g., 15 mL Wheaton); 8 volumes of solubilizing solution\*

\* The solubilizing solution is composed of 2% NP-40 (Sigma), 9 M urea (analytical grade, e.g., BDH or Bio-Rad), 0.5% dithiothreitol (DTT; Sigma) and 2% carrier ampholytes (pH 9-11 LKB: these come as a 20% stock solution, so 2% final concentration is achieved by making the final solution 10% 9-11 Ampholine by volume). A large batch of solubilizer (several hundred mL) is made and stored frozen at -80°C in aliquots sufficient to provide enough for one day's estimated sample preparation requirement. The solution is never allowed to become warmer than room temperature at any stage during preparation or thawing for use, since heating of concentrated urea solutions can produce contaminants that covalently modify proteins producing artifactual charge shifts. Once thawed, any unused solubilizer is discarded.

is added (i.e., 4 mL per 0.5 g tissue) and the mixture is homogenized using first the loose- and then the tight-fitting glass pestle. This takes approximately 5 strokes with each pestle and is carried out at room temperature because urea would crystallize out in the cold. Once the liver sample is thoroughly homogenized in the solubilizer, it is assumed that all the proteins are denatured (by the chaotropic effect of the urea and NP-40 detergent) and the enzymes inactivated by the high pH (~9.5). Therefore these samples may be kept at room temperature until they can be centrifuged or frozen as a group (within several hours of preparation). The samples are centrifuged for  $6 \times 10^4$  g min (e.g., 500 000  $\times$  g for 12 min using a Beckman TL-100 centrifuge). The centrifuge rotor is maintained at just below room temperature (e.g., 15–20°C), but not too cold, so as to prevent the precipitation of urea. The centrifuge of choice is a Beckman TL-100 because of the sample tube sizes available, but any ultracentrifuge accepting smallish tubes will suffice. When an appropriate centrifuge is not available near the site of sample preparation, samples can be frozen at –80°C and thawed prior to centrifugation and collection of supernatants. Each supernatant is carefully removed following centrifugation and aliquoted into at least 4 clean tubes for storage. This is done by transferring all the supernatant to one clean tube, mixing this gently (to assure homogeneous composition) and then dividing it into 4 aliquots. The aliquots are frozen immediately at –80°C. These multiple aliquots can provide insurance against a failed run or a freezer breakdown.

## 2.2 Two-dimensional electrophoresis

Sample proteins are resolved by 2-D electrophoresis using the 20  $\times$  25 cm Iso-Dalt® 2-D gel system ([26–29]; produced by LSB and by Hoefer Scientific Instruments, San Francisco) operating with 20 gels per batch. All first-dimensional isoelectric focusing (IEF) gels are prepared using the same single standardized batch of carrier ampholytes (BDH 4–8A in the present case, selected by LSB's batch-testing program for rat and mouse database work\*\*). A 10  $\mu$ L sample of solubilized liver protein is applied to each gel, and the gels are run for 33 000 to 34 500 volt-hours using a progressively increasing voltage protocol implemented by a programmable high-voltage power supply. An Angeliq™ computer-controlled gradient-casting system (produced by LSB) is used to prepare second-dimensional sodium dodecyl sulfate (SDS) polyacrylamide gradient slab gels in which the top 5% of the gel is 11%T acrylamide, and the lower 95% of the gel varies linearly from 11% to 18%T.

This system has recently been modified so as to employ a commercially available 30.8%T acrylamide/*N,N*'-methylenebisacrylamide prepared solution (thus avoiding the handling of the solid acrylamide monomer) and three additional stock solutions: buffer (made from Sigma pre-set Tris), persulfate and *N,N,N,N*'-tetramethylethylenediamine (TEMED). Each gel is identified by a computer-printed filter paper label polymerized into the lower left corner of the gel. First-dimensional IEF tube gels are loaded

directly (as extruded) onto the slab gels without equilibration, and held in place by polyester fabric wedges (Wedgies™, produced by LSB) to avoid the use of hot agarose. Second-dimensional slab gels are run overnight, in groups of 20, in cooled DALT tanks (10°C) with buffer circulation. All run parameters, reagent source and lot information, and notations of deviation from expected results are entered by the technician responsible on a detailed, multi-page record of the experiment.

## 2.3 Staining

Following SDS-electrophoresis, slab gels are stained for protein using a colloidal Coomassie Blue G-250 procedure in covered plastic boxes, with 10 gels (totalling approximately 1 L of gel) per box. This procedure (based on the work of Neuhauff [30, 31]) involves fixation in 1.5 L of 50% ethanol and 2% phosphoric acid for 2 h, three 30 min washes, each in 2 L of cold tap water, and transfer to 1.5 L of 34% methanol, 17% ammonium sulfate and 2% phosphoric acid for 1 h, followed by the addition of a gram of powdered Coomassie Blue G-250 stain. Staining requires approximately 4 days to reach equilibrium intensity, whereupon gels are transferred to cool tap water and their surfaces rinsed to remove any particulate stain prior to scanning. Gels may be kept for several months in water with added sodium azide. The water washes remove ethanol that would dissolve the stain (and render the system noncolloidal, with high backgrounds). The concentrated ammonium sulfate and methanol solution is diluted by equilibration with the water volume of the gels to automatically achieve the correct final concentrations for colloidal staining. Practical advantages of this staining approach can be summarized as follows: (i) the low, flat background makes computer evaluation of small spots (max OD < 0.02) possible, especially when using laser densitometry; (ii) up to 1500 spots can be reliably detected on many gels (e.g., rat liver) at loadings low enough to preserve excellent resolution; and (iii) reproducibility appears to be very good: at least several hundred spots have coefficients of reproducibility less than 15%. This value is at least as good as previous CBB methods, and significantly better than many silver stain systems.

## 2.4 Positional standardization

The carbamylated rabbit muscle creatine phosphokinase (CPK) standards [32] are purchased from Pharmacia and BDH. Amino acid compositions, and numbers of residues present in proteins used for internal standardization, are taken from the Protein Identification Resource (PIR) sequence database [33].

## 2.5 Computer analysis

Stained slab gels are digitized in red light at 134 micron resolution, using either a Molecular Dynamics laser scanner (with pixel sampling) or an Eikonix 78/99 CCD scanner. Raw digitized gel images are archived on high-density DAT tape (or equivalent storage media) and a greyscale video-print prepared from the raw digital image as hard-copy backup of the gel image. Gels are processed using the Kepler® software system (produced by LSB), a commercially available workstation-based software package built on

\*\* This material (succeeding certified batches of which are available from Hoefer Scientific Instruments) has the most linear pH gradient produced by any ampholyte tested except for the Pharmacia wide range (which has an unacceptable tendency to bind high-molecular weight acidic proteins, causing them to streak).

some of the principles of the earlier TYCHO system [34-41]. Procedure PROC008 is used to yield a spotlist giving position, shape and density information for each detected spot. This procedure makes use of digital filtering, mathematical morphology techniques and digital masking to remove the background, and uses full 2-D least-squares optimization to refine the parameters of a 2-D Gaussian shape for each spot. Processing parameters and file locations are stored in a relational database, while various log files detailing operation of the automatic analysis software are archived with the reduced data. The computed resolution and level of Gaussian convergence of each gel are inspected and archived for quality control purposes.

Experiment packages are constructed using the Kepler experiment definition database to assemble groups of 2-D patterns corresponding to the experimental groups (e.g., treated and control animals). Each 2-D pattern is matched to the appropriate "master" 2-D pattern (pattern F344MST3 in the case of Fischer 344 rat liver), thereby providing linkage to the existing rodent protein 2-D databases. The software allows experiments containing hundreds of gels to be constructed and analyzed as a unit, with up to 100 gels displayed on the screen at one time for comparative purposes and multiple pages to accommodate experiments of > 1000 gels. For each treatment, proteins showing significant quantitative differences vs. appropriate controls are selected using group-wise statistical parameters (e.g., Student's *t*-test, Kepler<sup>®</sup> procedure STUDENT). Proteins satisfying various quantitative criteria (such as  $P < 0.001$  difference from appropriate controls) are represented as highlighted spots onscreen or on computer-plotted protein maps and stored as spot populations (*i.e.*, logical vectors) in a liver protein database. Quantitative data (spot parameters, statistical or other computed values) are stored as real-valued vectors in the database. Analysis of coregulation is performed using a Pierson product-moment correlation (Kepler procedure CORREL) to determine whether groups of proteins are coordinately regulated by any of the treatments. Such groups can be presented graphically on a protein map, and reported together with the statistical criteria used to assess the level of coregulation. Multivariate statistical analysis (e.g., principal components' analysis) is performed on data exported to SAS (SAS Institute).

## 2.6 Graphical data output

Graphical results are prepared in GKS and translated within Kepler<sup>®</sup> into output for any of a variety of devices. Linedrawing output is typically prepared as Postscript and printed on an Apple Laserwriter. Detailed maps presented here have been generated using an ultra-high-resolution Postscript-compatible Linotronic output device. Greyscale graphics are reproduced from the workstation screen using a Seikosha videoprinter. Patterns are shown in the standard orientation, with high molecular mass at the top and acidic proteins to the left.

## 2.7 Experiment LSBC04

In the study described here 12-week-old Charles River male F344 rats were used. Diets were prepared at LSB, based on a Purina 5755M Basal Purified Diet. Lovastatin and cholestyramine were obtained as prescription pharma-

ceuticals, ground and mixed with the diet at concentrations of 0.075% and 1%, respectively. The high cholesterol diet was Purina 5801M-A (5% cholesterol plus 1% sodium cholate in the control diet). Animal work was carried out by Microbiological Associates (Bethesda, MD). Animals were acclimatized for one week on the control diet, fed test or control diets for one week, and sacrificed on day 8. Average daily doses of lovastatin and cholestyramine in appropriate groups were 37 mg/kg/day and 5 g/kg/day, respectively, based on the weight of the food consumed. Liver samples were collected and prepared for 2-D electrophoresis according to the standard liver protocol (homogenization in 8 volumes of 9 M urea, 2% NP-40, 0.5% dithiothreitol, 2% LKB pH 9-11 carrier ampholytes, followed by centrifugation for 30 min at 80 000  $\times$  g). Kidney, brain and plasma samples were frozen. Gels were run as described above, and the data was analyzed using the Kepler<sup>®</sup> system. Gels were scaled, to remove the effect of differences in protein loading, by setting the summed abundances of a large number of matched spots equal for each gel (linear scaling).

## 3 Results and discussion

### 3.1 The rat liver protein 2-D map

F344MST3 is a standard 2-D pattern of rat liver proteins, based on the Fischer 344 strain. This pattern was initiated from a single 2-D gel and extensively edited in an experiment comparing it to a range of protein loads, so as to include both small spots and well-resolved representations of high-abundance spots. More than 700 rat liver 2-D patterns have been matched to F344MST3 in a series of drug effects and protein characterization experiments, and numerous new spots (induced by specific drugs, for instance) have been added as a result. A modified version including additional spots present in the Sprague-Dawley outbred rat has also been developed (data not shown). Figure 1 shows a greyscale representation and Fig. 2 a schematic plot of the master pattern. More than 1200 spots are included, most of which are visible on typical gels loaded with 10  $\mu$ L of solubilized liver protein prepared by the standard method and stained with colloidal Coomassie Blue. Master spot numbers (MSN's) have been assigned to all proteins, and appear in the following figures, each showing one quadrant of the pattern. Figure 3 shows the upper left (acidic, high molecular mass) quadrant, Fig. 4 the upper right (basic, high molecular mass) quadrant, Fig. 5 the lower left (acidic, low molecular mass) quadrant, and Fig. 6 the lower right (basic, low molecular mass) quadrant. The quadrants overlap as an aid to moving between them. The gel position (in 100 micron units), isoelectric point (relative to the CPK internal *p*/standards) and SDS molecular mass (from the calibration curve in Fig. 8) are listed for each spot (Table 1). Because of the precision of the CPK-*p*/values, these parameters can be used to relate spot locations between gel systems more reliably than using *p*/measurements expressed as pH. A major objective of current studies is the identification of all major spots corresponding to known liver proteins, as well as rigorous definitions of subcellular organelle contents. Of particular interest to us is the parallel development of identifications in the rat and mouse liver maps, allowing detailed comparisons of gene expression effects in the two systems. The results of these studies will be presented systematically in a later edition of this database,

but we include here a useful series of 22 orienting identifications as an aid to other users of the rat liver pattern (Table 2).

### 3.2 Carbamylated charge standards, computed *pI*'s and molecular mass standardization

We have previously shown that the use of a system of closely-spaced internal *pI* markers (made by carbamylating a basic protein) offers an accurate and workable solution to the problem of assigning positions in the *pI* dimension [32]. The same system, based on 36 protein species made by carbamylating rabbit muscle CPK, has been used here to assign *pI*'s to most rat liver acidic and neutral proteins. The standards were coelectrophoresed with total liver proteins, and the standard spots added to a special version of the master pattern F344MST3. The gel *X*-coordinates of all liver protein spots lying within the CPK charge train were then transformed into CPK *pI* positions by interpolation between the positions of immediately adjacent standards (Table 1) using a Kepler® vector procedure.

It has proven possible to compute fairly accurate *pI* values for many proteins from the amino acid composition [42]. We have attempted here to test a further elaboration of this approach, in which we computed *pI*'s for the CPK standards themselves, based on our knowledge of the rabbit muscle CPK sequence and the fact that adjacent members of the charge train typically differ by blockage of one additional lysine residue (Table 3). We compared these values to similar computed *pI*'s for an additional set of carbamylated standards made from human hemoglobin beta chains and a series of rat liver and human plasma proteins of known position and sequence (Fig. 7, Table 4). The result demonstrates good concordance between these systems. Two proteins show significant deviations: liver fatty-acid binding protein (FABP; #1 in Table 4) and protein disulphide isomerase (#20 in the table). The FABP spot present on F344MST3 may represent a charge-modified version of a more basic parent spot closer to the expected *pI*, not resolved in the IEF/SDS gel. Of particular importance is the fact that, by comparing computed *pI*'s of sequenced but unlocated proteins with the CPK *pI*'s, we can assign a probable gel location without making any assumptions regarding the actual gel pH gradient. This offers a useful shortcut, given the vagaries of pH measurement on small diameter IEF gels. We have used this approach to compute the CPK *pI*'s of all rat and mouse proteins in the PIR sequence database, as an aid to protein identification (data not shown).

In order to standardize SDS molecular weight (SDS-MW), we have used a standard curve fitted to a series of identified proteins (Fig. 8). Rather than using molecular mass *per se*, we have elected to use the number of amino acids in the polypeptide chain, as perhaps a better indication of the length of the SDS-coated rod that is sieved by the second dimension slab. The resulting values were multiplied by 112 (the weighted average mass of amino acids in sequenced proteins) to give predicted molecular masses. Because we use gradient slabs, we have not constrained the fitted curve to conform to any predetermined model; rather we tried many equations and selected the best using the program "Tablecurve" on a PC. The equation chosen was  $y = a + bx + c/x^2$ , where *y* is the number of residues, *x* is the gel

*Y*-coordinate, *a* is 511.83, *b* is -0.2731 and *c* is 33183801. The resulting fit appears to be fairly good over a broad range of molecular mass.

### 3.3 An example of rat liver gene regulation: Cholesterol metabolism

Experiment LSBC04 was designed as a small-scale test of the regulation of cholesterol metabolism *in vivo* by three agents included in the diet: lovastatin (Mevacor®, an inhibitor of HMG-CoA reductase); cholestyramine (a bile acid sequestrant that has the effect of removing cholesterol from the gut-liver recirculation); and cholesterol itself. The first two agents should lower available cholesterol and the third should raise it, allowing manipulation of relevant gene expression control systems in both directions. Such an experiment offers an interesting test of the 2-D mapping system since most of the pathway enzymes are present in low abundance, many are membrane-bound and difficult to solubilize, and the pathway itself is complex. Approximately 1000 proteins were separated and detected in liver homogenates. Twenty-one proteins were found to be affected by at least one treatment, and these could be divided into several coregulated groups.

#### 3.3.1 MSN 413 (putative cytosolic HMG-CoA synthase) and sets of spots regulated coordinately or inversely

One group of spots (including a spot assigned to the cytosolic HMG-CoA synthase, MSN 413) showed the expected increase in abundance with lovastatin or cholestyramine, the synergistic further increase with lovastatin and cholestyramine, and a dramatic decrease with the high cholesterol diet. Spot number 413 is the most strongly regulated protein in the present experiment, showing a 5- to 10-fold induction after a 1 week treatment with 0.075% lovastatin and 1% cholestyramine in the diet (Figs. 9 and 10). Its expression follows precisely the expectation for an enzyme whose abundance is controlled by the cholesterol level; it is progressively increased from the control levels by cholestyramine, lovastatin and lovastatin plus cholestyramine, and it sinks below the threshold of detection in animals fed the high cholesterol diet. This spot has been tentatively identified as the cytosolic HMG-CoA synthase, based on a reaction with an antiserum to that protein provided by Dr. Michael Greenspan at Merck Sharp & Dohme Research Laboratories. This enzyme lies immediately before HMG-CoA reductase in the liver cholesterol biosynthesis pathway, and is known to be co-regulated with it. Spot 413 has an SDS molecular weight of about 54 000 and a CPK *pI* of -11.4, in reasonably close agreement with a molecular weight of 57 300 and a CPK *pI* of -15.7 computed from the known sequence of the hamster enzyme [43].

Using a classical product-moment correlation test (Kepler procedure CORREL), a series of five additional spots was found to be coregulated with 413. The level of correlation was exceedingly high (> 95%). Two of these, 1250 and 933, are at similar molecular weights and approximately one charge more acidic than 413 (Fig. 9), indicating that they may be covalently modified forms of the 413 polypeptide. This suspicion is strengthened by the observation that both spots are also stained by the antibody to cytosolic HMG-CoA synthase. The remaining three correlated spots appear

to comprise an additional related pair (1253 and 1001) of around 40 kDa and a single spot (1119) of around 28 kDa. Because these two presumed proteins are present at substantially lower abundances than 413, and because the cytosolic HMG-CoA synthase is reported to consist of only one type of polypeptide, they are likely to represent other, very tightly coregulated enzymes. A second group of six spots was selected based on a regulatory pattern close to the inverse of that for spot 413 (MSN's 34, 79, 178, 182, 204, 347; data not shown). For these proteins, the lowest level of expression occurs with exposure to lovastatin plus cholestyramine and the highest level upon exposure to the high-cholesterol diet. Spots 182 and 79 are highly correlated and lie about one charge apart at the same molecular weight; they may thus be isoforms of a single protein. The other four spots probably represent additional enzymes or subunits.

### 3.3.2 MSN 235 and coregulated spots

A third group of five spots, mainly comprised of mitochondrial proteins including putative mitochondrial HMG-CoA synthase spots, showed a modest induction by lovastatin alone, but little or no effect with any of the other treatments (including the combination of lovastatin and cholestyramine; Fig. 12). This result is intriguing because lovastatin was expected to affect only the regulation of enzymes of cholesterol synthesis, which is entirely extra-mitochondrial. Three of the spots (235, 134, 144) form a closely-packed triad at approximately 30 kDa, and are likely to represent isoforms of one protein. All three spots are stained by an antibody to the mitochondrial form of HMG-CoA synthase obtained from Dr. Greenspan. Subcellular fractionation indicates a mitochondrial location. The other two spots (633 at about 38 kDa and 724 at about 69 kDa) are each present at lower abundance than the members of the triad.

### 3.3.3 An example of an anti-synergistic effect

A sixth spot (367) shows strong induction by lovastatin (two- to threefold), and about half as much induction with lovastatin plus cholestyramine, but without sharing the animal-animal heterogeneity pattern of the 235-set (Fig. 13). This protein is also mitochondrial, and represents the clearest example of an anti-synergistic effect of lovastatin and cholestyramine. The existence of such an effect demonstrates that lovastatin and cholestyramine do not act exclusively through the same regulatory pathway.

### 3.3.4 Complexity of the cholesterol synthesis pathway

Taken together, these results suggest that treatment with lovastatin alone can affect both cytosolic and mitochondrial pathways using HMG-CoA, while cholestyramine, on the other hand, either alone or in combination with lovastatin, produces a strong effect on the putative cytosolic pathway, but little or no effect on the putative mitochondrial pathway. An explanation for this difference may lie in lovastatin's effect on levels of HMG-CoA and related precursor compounds that are exchanged between the cytosol and the mitochondrion, whereas cholestyramine should affect only the cytosolic pathways directly controlled by cholesterol and bile acid levels. It remains to be explained why some

proteins of the putative mitochondrial pathway are so much more variable in their expression in all groups. An examination of all the coregulated groups suggests that quantitative statistical techniques can extract a wealth of interesting information from large sets of reproducible gels. The abundance of spots in the 413 coregulation group, for example, shows an amazing level of concordance in their relative expression among the five individuals of the lovastatin and cholestyramine treatment group. This effect is not due to differences in total protein loading, since they have already been removed by scaling, and since proteins with quite different regulation patterns can be demonstrated (e.g., Fig. 13). Such effects raise the possibility that many gene coregulation sets may be revealed through the study of a sufficiently large population of control animals (i.e., without any experimental manipulation). This approach, exploiting natural biological variation in protein expression instead of drug effects, offers an important incentive for the construction of a large library of control animal patterns.

## 4 Conclusions

Because of the widespread use of rat liver in both basic biochemistry and in toxicology, there is a long-term need for a comprehensive database of liver proteins. The rat liver master pattern presented here has proven to be an accurate representation of this system, having been matched to more than 700 gels to date. As the number of proteins identified and the number of compounds tested for gene expression effects grows, we expect this database to contribute valuable insights into gene regulation. Its practical utility in several areas of mechanistic toxicology is already being demonstrated.

Received September 11, 1991

## 5 References

- [1] O'Farrell, P., *J. Biol. Chem.* 1975, 250, 4007-4021.
- [2] Klose, J., *Humangenetik* 1975, 26, 231-243.
- [3] Scheele, G. A., *J. Biol. Chem.* 1975, 250, 5375-5385.
- [4] Iborra, G. and Buhler, J. M., *Anal. Biochem.* 1976, 74, 503-511.
- [5] Anderson, N. G. and Anderson, N. L., *Behring. Inst. Mitt.* 1979, 63, 169-210.
- [6] Anderson, N. G. and Anderson, N. L., *Clin. Chem.* 1982, 28, 739-748.
- [7] Heydorn, W. E., Creed, G. J. and Jacobowitz, D. M., *J. Pharmacol. Exp. Therap.* 1984, 229, 622-628.
- [8] Anderson, N. L., Nance, S. L., Tollaksen, S. L., Giere, F. A. and Anderson, N. G., *Electrophoresis* 1985, 6, 592-599.
- [9] Racine, R. R. and Langley, C. H., *Biochem. Genet.* 1980, 18, 185-197.
- [10] Klose, J., *Mol. Evol.* 1982, 18, 315-328.
- [11] Neel, J. V., Baier, L., Hanash, S. and Erickson, R. P., *J. Hered.* 1985, 76, 314-320.
- [12] Marshall, R. R., Raj, A. S., Grant, F. J. and Heddle, J. A., *Can. J. Genet. Cytol.* 1983, 25, 457-446.
- [13] Taylor, J., Anderson, N. L., Anderson, N. G., Gemmell, A., Giometti, C. S., Nance, S. L. and Tollaksen, S. L., in: Dunn, M. J. (Ed.), *Electrophoresis '86*, Verlag Chemie, Weinheim 1986, pp. 583-587.
- [14] Giometti, C. S., Gemmell, M. A., Nance, S. L., Tollaksen, S. L. and Taylor, J., *J. Biol. Chem.* 1987, 262, 12764-12767.
- [15] Anderson, N. L., Giere, F. A., Nance, S. L., Gemmell, M. A., Tollaksen, S. L. and Anderson, N. G., in: Galteau, M.-M. and Siest, G. (Eds.), *Progrès Récents en Electrophorèse Bidimensionnelle*, Presses Universitaires de Nancy, Nancy 1986, pp. 253-260.
- [16] Anderson, N. L., Swanson, M., Giere, F. A., Tollaksen, S., Gemmell, A., Nance, S. L. and Anderson, N. G., *Electrophoresis* 1986, 7, 44-48.

- [17] Anderson, N. L., Giere, F. A., Nance, S. L., Gemmell, M. A., Tollaksen, S. L. and Anderson, N. G., *Fundam. Appl. Toxicol.* 1987, 8, 39-50.
- [18] Anderson, N. L., in: *New Horizons in Toxicology*, Eli Lilly Symposium, 1991, in press.
- [19] Antoine, B., Rahimi-Pour, A., Siest, G., Magdalou, J. and Galteau, M. M., *Cell. Biochem. Funct.* 1987, 5, 217-231.
- [20] Elliott, B. M., Ramasamy, R., Stonard, M. D. and Spragg, S. P. *Biochim. Biophys. Acta* 1986, 870, 135-140.
- [21] Huber, B. E., Heilman, C. A., Wirth, P. J., Miller, M. J. and Thorpeirson, S. S., *Hepatology* 1986, 6, 209-219.
- [22] Wirth, P. J. and Vesterberg, O., *Electrophoresis* 1988, 9, 47-53.
- [23] Witzmann, F. A. and Parker, D. N., *Toxicol. Lett.* 1991, 57, 29-36.
- [24] Rampersaud, A., Waxman, D. J., Ryan, D. E., Levin, W. and Walz, F. G., Jr., *Arch. Biochem. Biophys.* 1985, 243, 174-183.
- [25] Vlasuk, G. P. and Walz, F. G., Jr., *Anal. Biochem.* 1980, 105, 112-120.
- [26] Anderson, N. G. and Anderson, N. L., *Anal. Biochem.* 1978, 85, 331-340.
- [27] Anderson, N. L. and Anderson, N. G., *Anal. Biochem.* 1978, 85, 341-354.
- [28] Anderson, L., Hofmann, J.-P., Anderson, E., Walker, B. and Anderson, N. G., in: Endler, A. T. and Hanash, S. (Eds.), *Two-Dimensional Electrophoresis*, VCH Verlagsgesellschaft, Weinheim 1989, pp. 288-297.
- [29] Anderson, L., *Two-Dimensional Electrophoresis: Operation of the ISO-DALT<sup>®</sup> System*, Large Scale Biology Press, Washington, DC 1988, ISBN 0-945532-00-8, 170pp.
- [30] Neuhoft, V., Stamm, R. and Eibl, H., *Electrophoresis* 1985, 6, 427-448.
- [31] Neuhoft, V., Arold, N., Taube, D. and Ehrhardt, W., *Electrophoresis* 1988, 9, 255-262.
- [32] Anderson, N. L. and Hickman, B. J., *Anal. Biochem.* 1979, 93, 312-320.
- [33] Sidman, K. E., George, D. E., Barker, W. C. and Hunt, L. T., *Nucl. Acids Res.* 1988, 16, 1869-1871.
- [34] Taylor, J., Anderson, N. L., Coulter, B. P., Scandora, A. E. and Anderson, N. G., in: Radola, B. J. (Ed.), *Electrophoresis '79*, de Gruyter, Berlin 1980, pp. 329-339.
- [35] Taylor, J., Anderson, N. L. and Anderson, N. G., in: Allen, R. C. and Arnaud, P. (Eds.), *Electrophoresis '81*, de Gruyter, Berlin 1981, pp. 383-400.
- [36] Anderson, N. L., Taylor, J., Scandora, A. E., Coulter, B. P. and Anderson, N. G., *Clin. Chem.* 1981, 27, 1807-1820.
- [37] Taylor, J., Anderson, N. L., Scandora, A. E., Jr., Willard, K. E. and Anderson, N. G., *Clin. Chem.* 1982, 28, 861-866.
- [38] Taylor, J., Anderson, N. L. and Anderson, N. G., *Electrophoresis* 1983, 4, 338-345.
- [39] Anderson, N. L. and Taylor, J., in: *Proceedings of the Fourth Annual Conference and Exposition of the National Computer Graphics Association*, Chicago, June 26-30, 1983, pp. 69-76.
- [40] Anderson, N. L., Hofmann, J.-P., Gemmell, A. and Taylor, J., *Clin. Chem.* 1984, 30, 2031-2036.
- [41] Anderson, L., in: Schafer-Nielsen, C. (Ed.), *Electrophoresis '88*, VCH Verlagsgesellschaft, Weinheim 1988, pp. 313-321.
- [42] Neidhardt, F. C., Appleby, D. A., Sankar, P., Hutton, M. E. and Phillips, T. A., *Electrophoresis* 1989, 10, 116-121.
- [43] Gil, G., Goldstein, J. L., Slaughter, C. A. and Brown, M. S., *J. Biol. Chem.* 1986, 261, 3710-3716.

## 6 Addendum 1: Figures 1-13

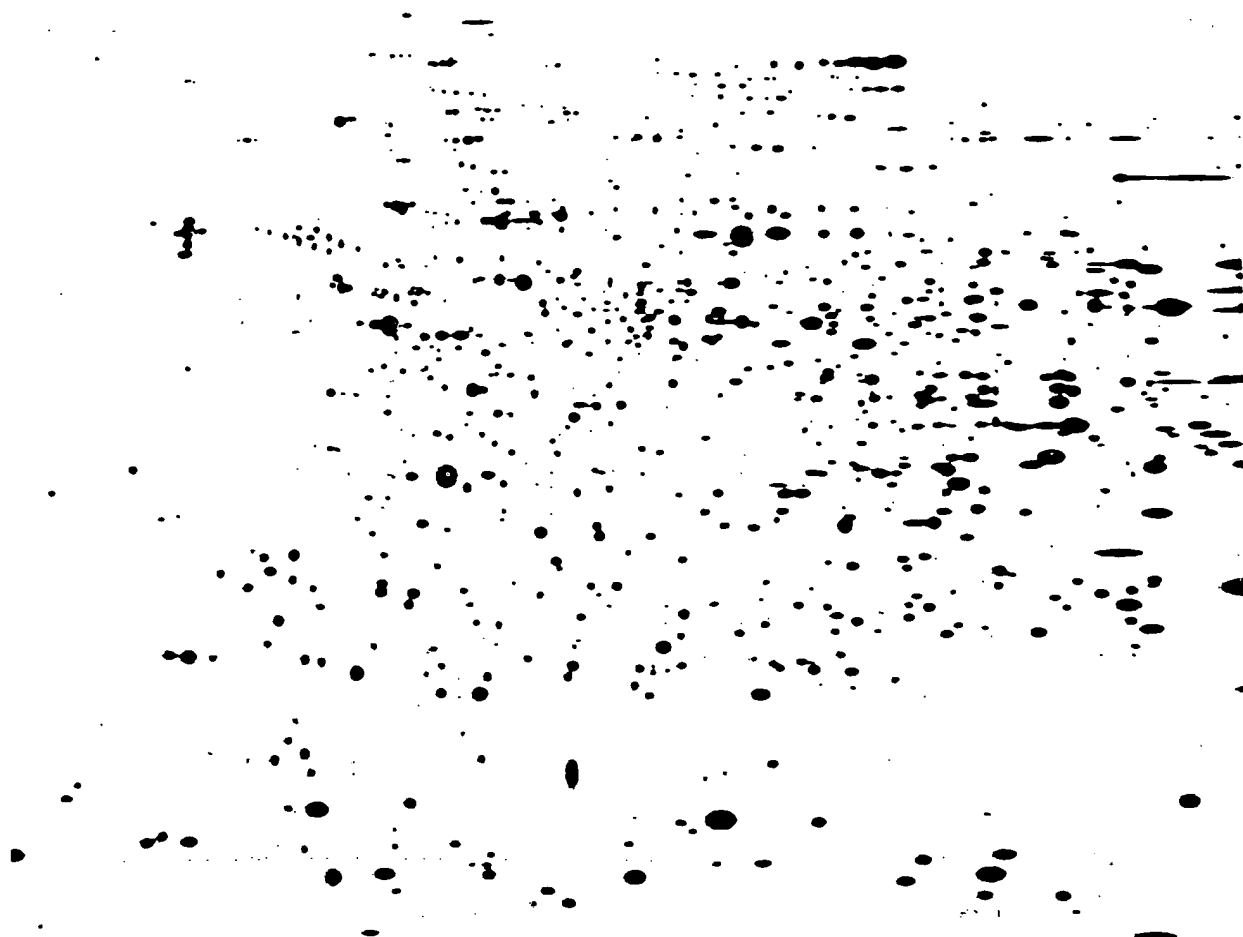


Figure 1. Synthetic representation of the standard rat liver 2-D master pattern, rendered as a greyscale image using a videoprinter.

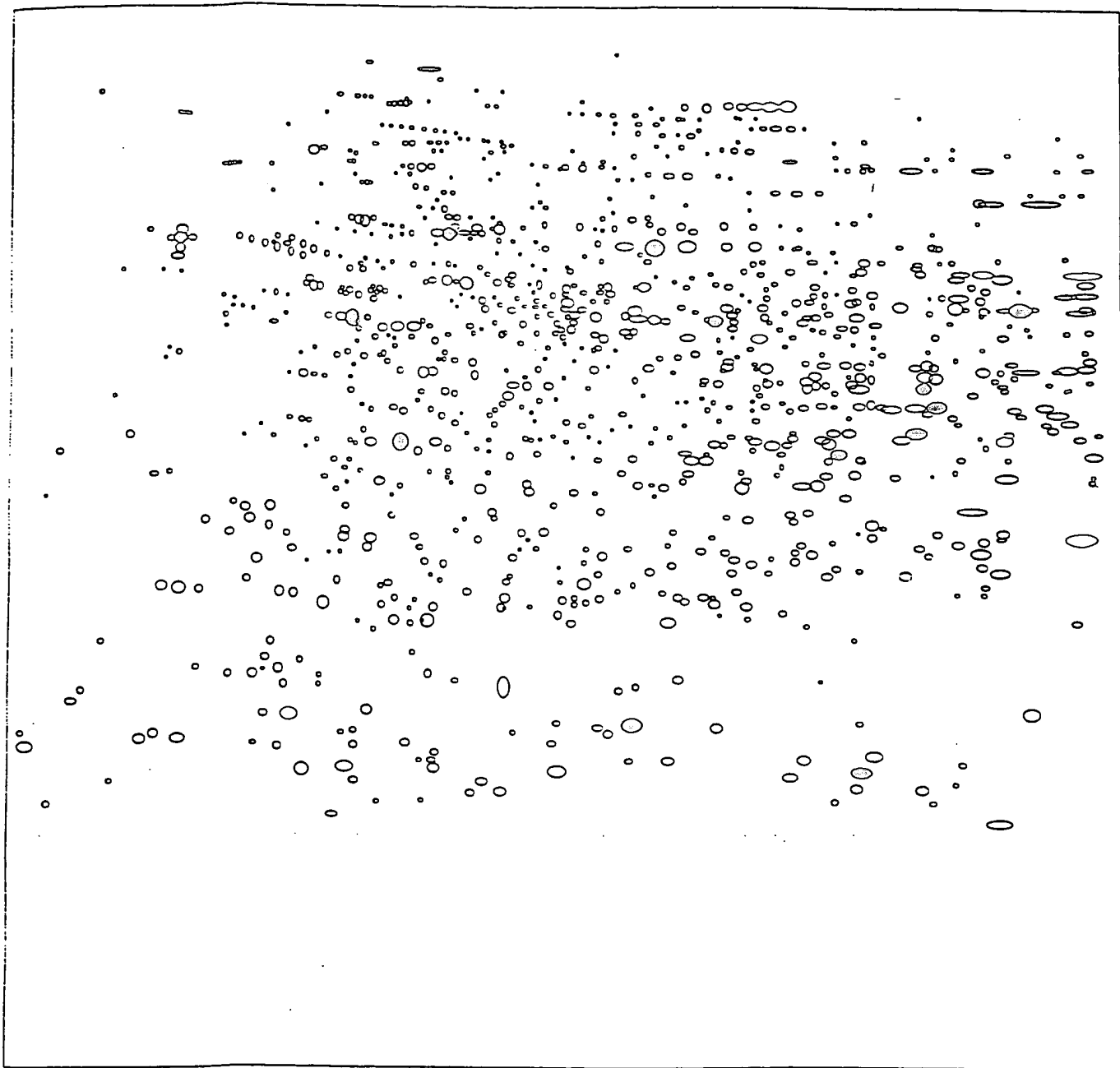


Figure 2. Schematic representation of the master pattern (the same as Fig. 1), useful as an aid in relating specific areas of Fig. 1 and the following detailed quadrants.

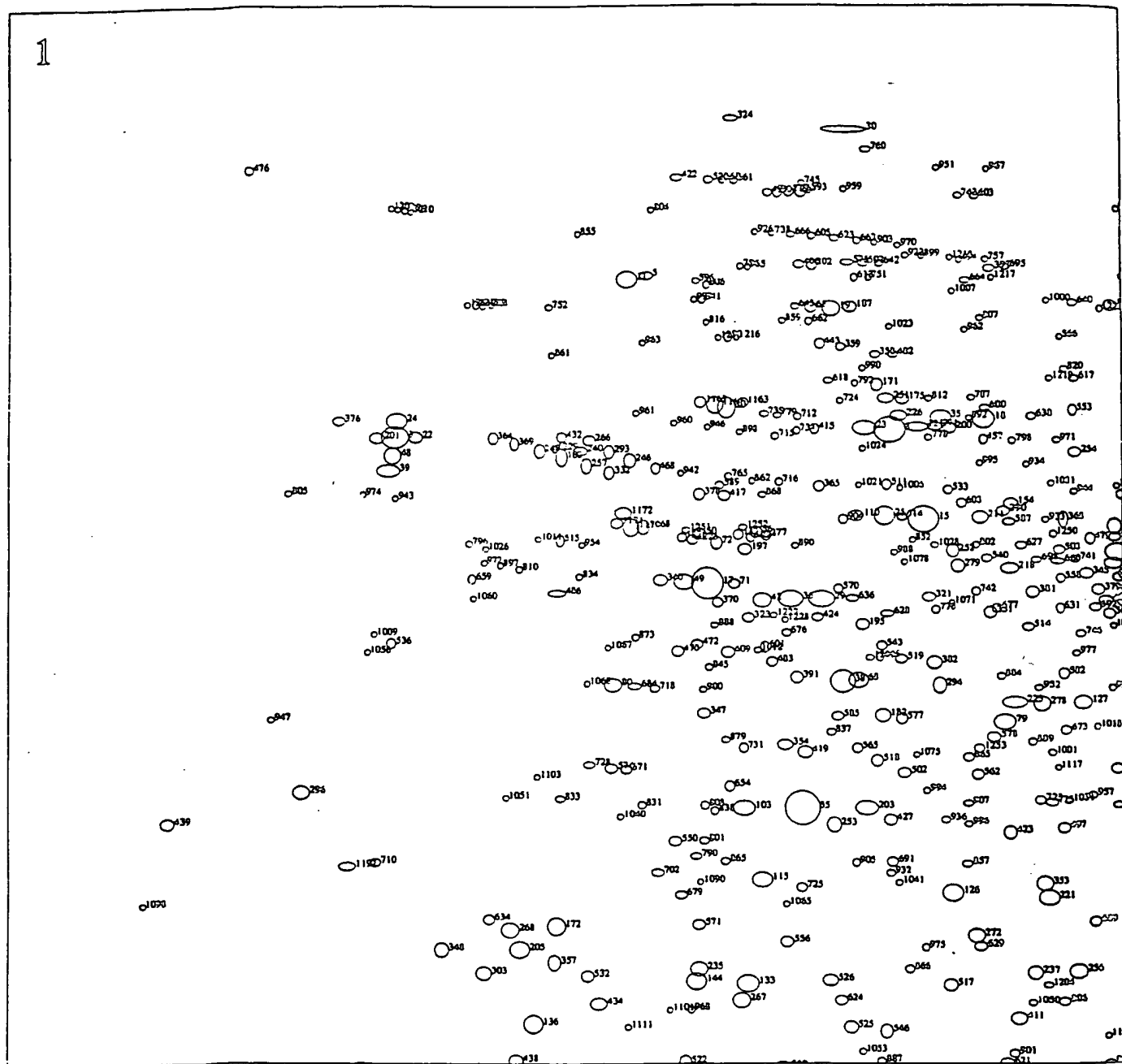


Figure 3. Upper left (high molecular weight, acidic) quadrant (#1) of the rat liver map, showing spot numbers.

2

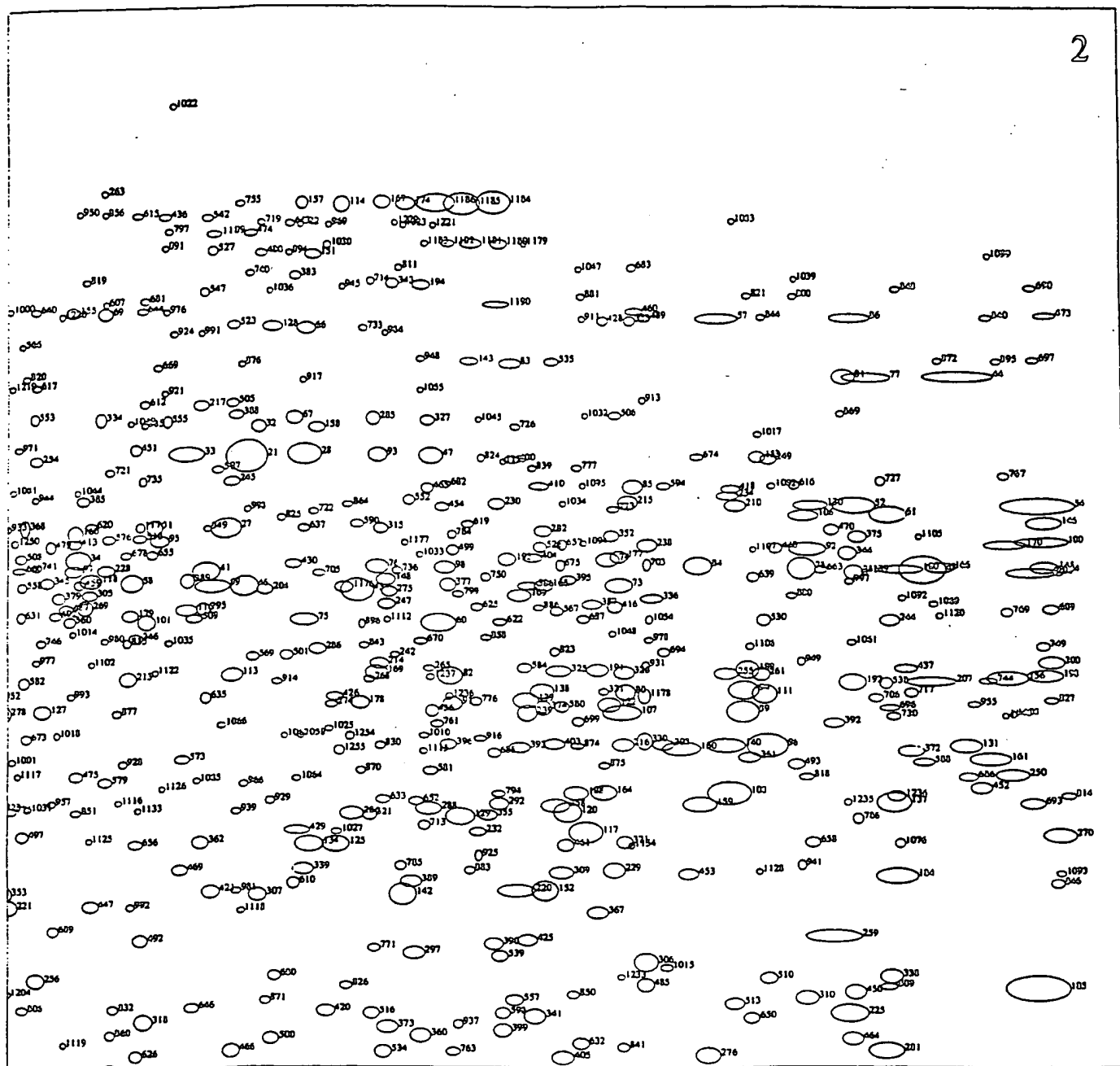


Figure 4. Upper right (high molecular weight, basic) quadrant (#2) of the rat liver map, showing spot numbers.

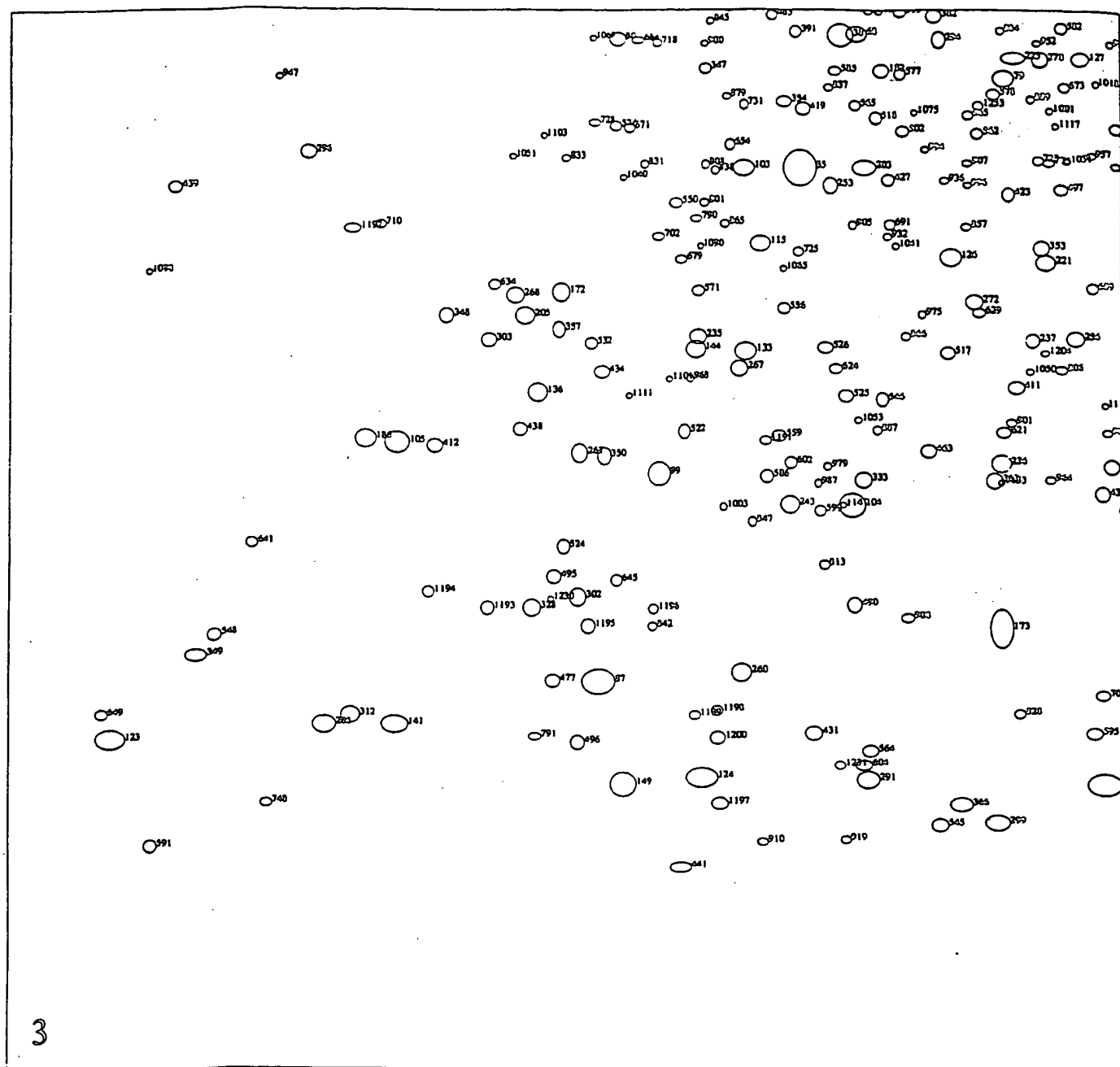


Figure 5. Lower left (low molecular weight, acidic) quadrant (#3) of the rat liver map, showing spot numbers.

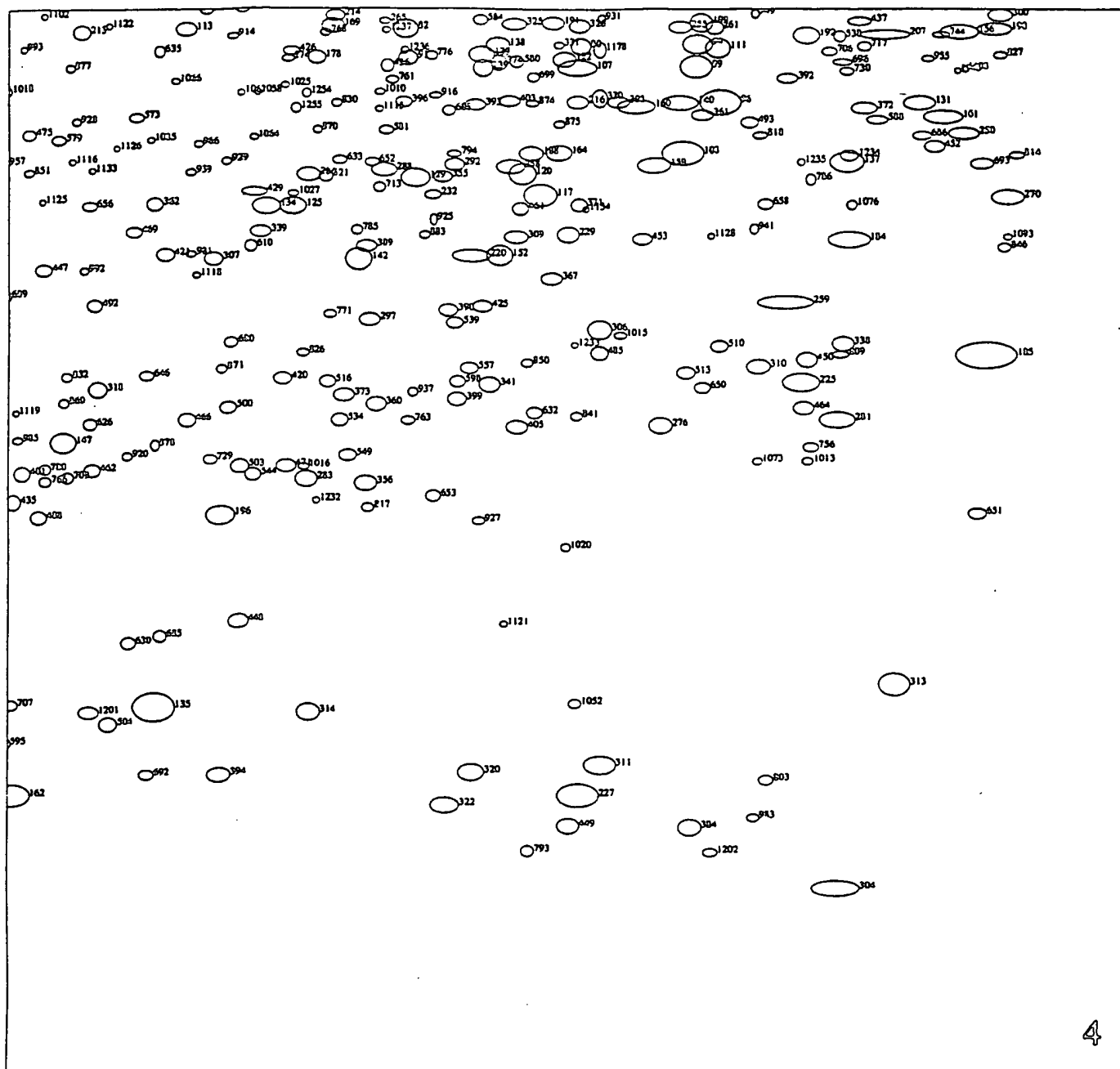


Figure 6. Lower right (low molecular weight, basic) quadrant (#4) of the rat liver map, showing spot numbers.

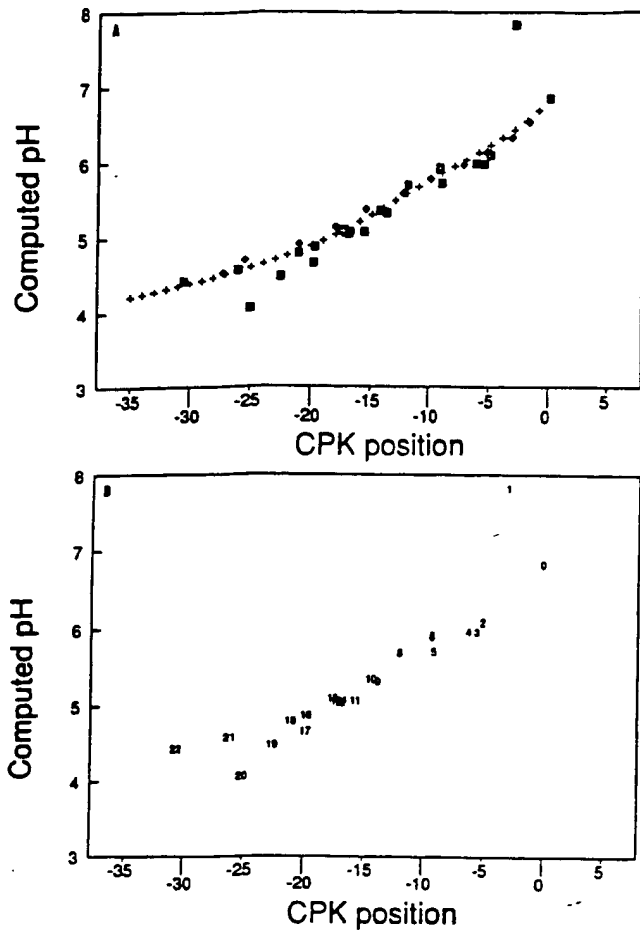


Figure 7. (a) Plot of computed isoelectric point versus gel X-position for two sets of carbamylated standard proteins (rabbit muscle CPK [+]) and human hemoglobin  $\beta$  chain, filled diamonds) and several other proteins (shaded squares). (b) The identities of the various proteins represented by the squares are indicated by the numbers in corresponding positions on (a); these refer to Table 4.

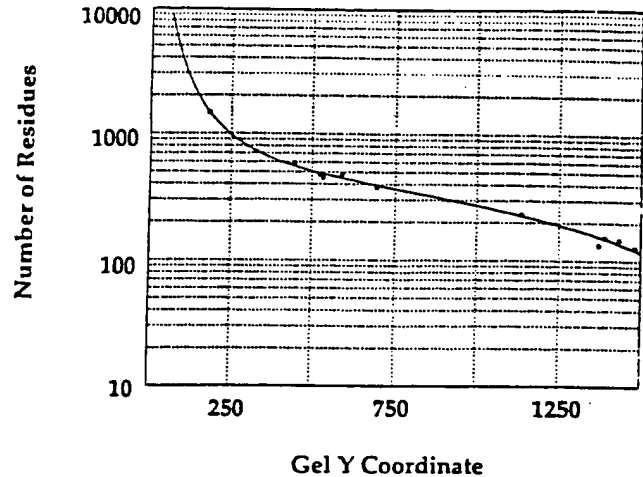


Figure 8. Plot of number of amino acids versus gel Y-position, with fitted curve used to predict molecular mass of unidentified proteins.

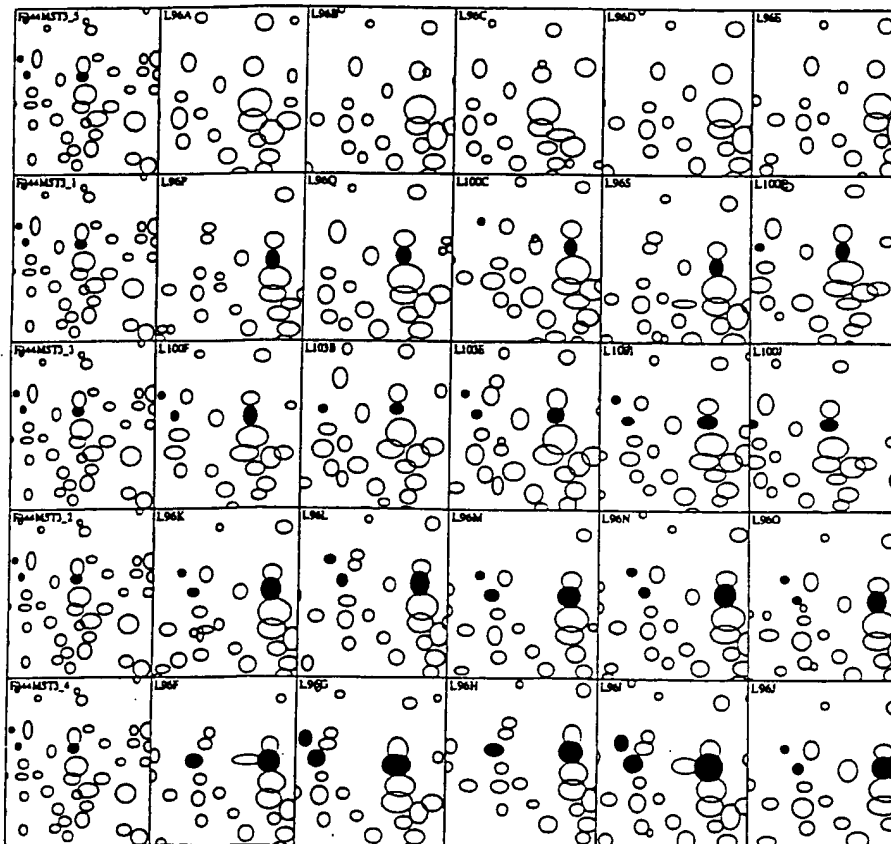


Figure 9. Montage showing effects in the region of MSN:413. The montage shows a small window into one portion of the 2-D pattern, one row of windows for each experimental group, and one panel for each gel in the experiment. The left-most pattern in each row is a group-specific copy of the master pattern followed by the patterns for the five individual rats in the group. The highlighted protein spots (filled circles) are spot 413 (on the right of each panel; identified as cytosolic HMG-CoA synthase) and two modified forms of it (1250 and 933). From the top, the rows (experimental groups) are: high cholesterol, controls, cholestyramine, lovastatin, and lovastatin plus cholestyramine.

## Regulation of Rat Liver 413

(Putative Cytosolic HMG-CoA Synthase, 53kd)  
Test Compounds in Diet

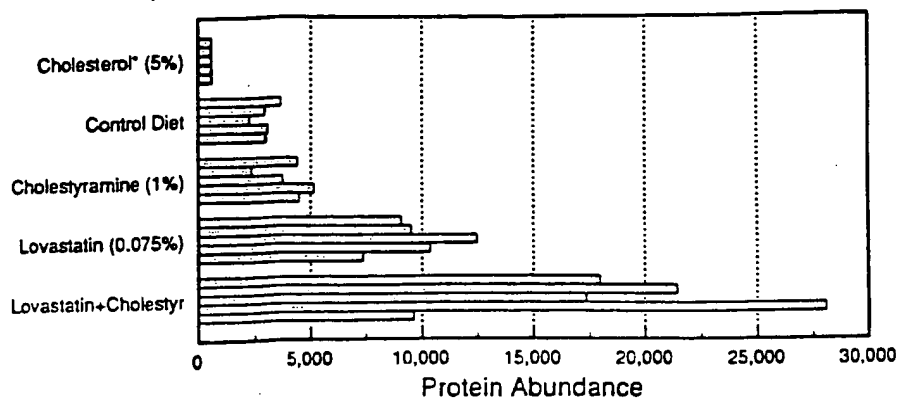


Figure 10. Bargraph showing the quantitative effects of various treatments on the abundance of MSN:413 (cytosolic HMG-CoA synthase) in the gels of Fig. 9.

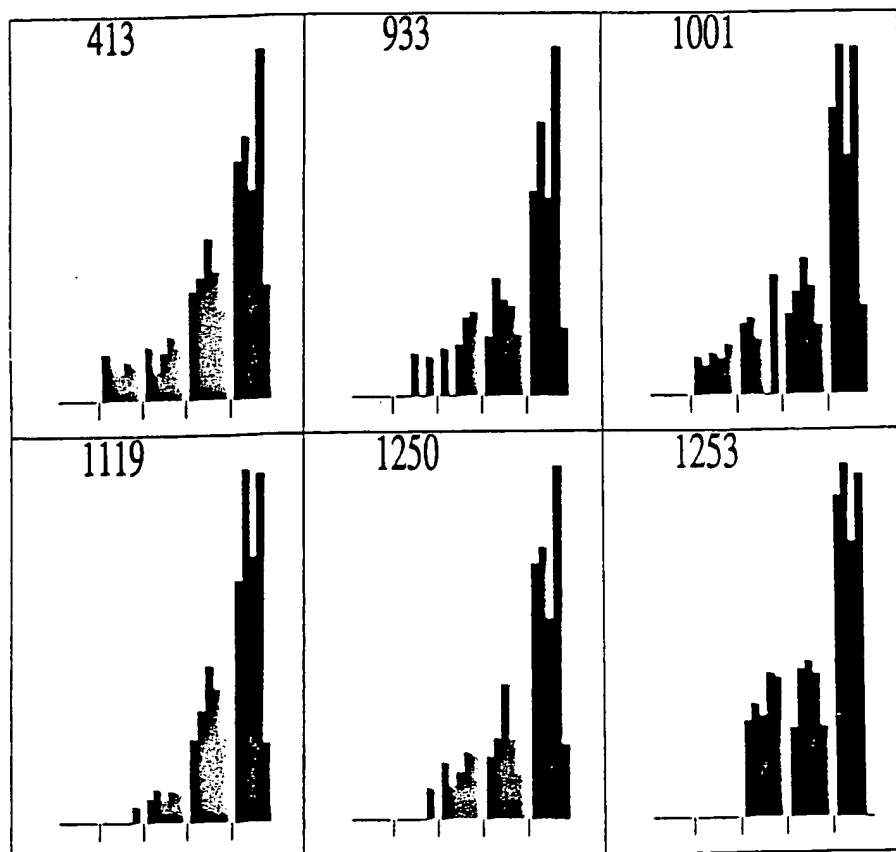


Figure 11. Bargraphs of a series of six coregulated spots including MSN:413. In the bargraphs, the abundances of the appropriate spot (master spot number shown at the top of the panel) in each animal are shown. The five five-animal groups are in the order (left to right): high cholesterol, controls, cholestyramine, lovastatin, and lovastatin plus cholestyramine. Each bar within a group represents one experimental animal liver (one 2-D gel). Note the correlated expression of the 6 spots, especially in the two far right (most strongly induced) groups.

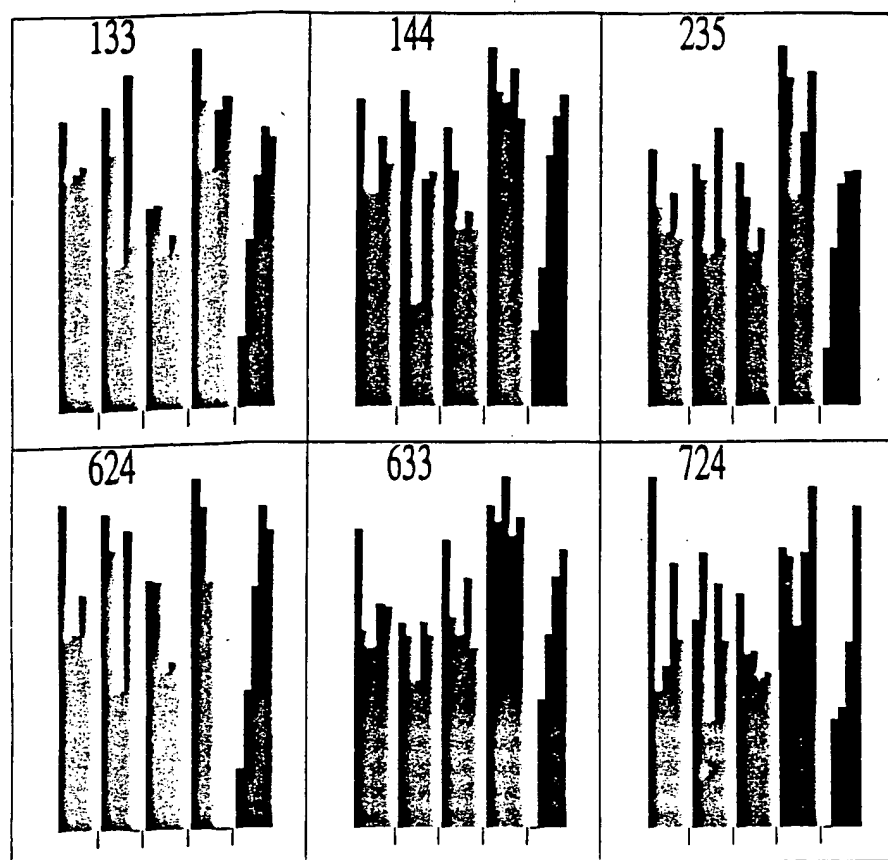


Figure 12. Data on a second coregulated group of spots, presented as in Fig. 11. The fourth experimental group (lovastatin) shows a modest induction, while the fifth group (lovastatin plus cholestyramine) does not.

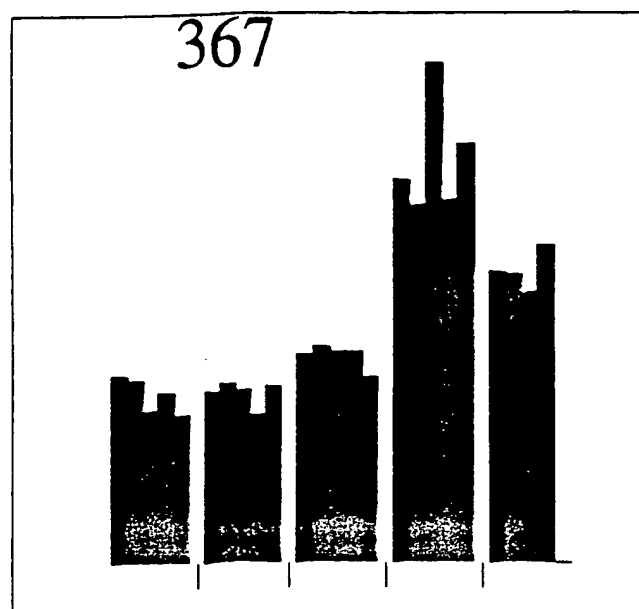


Figure 13. Data on spot MSN:367, presented as in Fig. 11. This protein shows unambiguously the anti-synergistic effect of lovastatin and cholestyramine (fifth group) as compared to lovastatin (fourth group). This response contrasts strongly with the regulation pattern seen in Fig. 11.

Table 1. Master table of proteins in the rat liver database<sup>a)</sup>

MSN	X	Y	CPKpI	SDSMW	MSN	X	Y	CPKpI	SDSMW	MSN	X	Y	CPKpI	SDSMW
3	311	434	<-35.0	63,800	95	1119	536	-9.9	53,800	174	1364	183	-6.7	162,900
5	568	263	-24.3	102,900	96	1731	756	-2.0	40,700	175	825	393	-15.7	69,300
8	812	426	-16.0	64,800	97	1033	566	-11.4	51,600	177	1582	553	-3.6	52,600
11	549	268	-25.2	101,000	98	1406	565	-6.1	51,700	178	1321	710	-7.2	43,000
15	845	520	-15.3	55,200	99	578	1149	-23.8	25,000	179	1089	615	-10.4	48,300
17	629	589	-21.6	50,000	100	2004	538	>0.0	53,700	180	1866	567	-0.5	51,600
18	906	414	-14.0	66,300	101	1106	623	-10.1	47,900	181	411	295	-32.1	91,200
19	755	298	-17.5	90,200	102	482	455	-28.5	61,300	182	804	730	-16.2	42,000
20	649	403	-20.9	67,900	103	665	830	-20.2	37,300	184	1860	896	-0.6	34,500
21	1204	448	-8.7	62,100	104	773	1182	-17.0	23,800	185	1997	1017	>0.0	29,800
22	332	434	<-35.0	63,800	105	312	1117	<-35.0	26,100	186	279	1113	<-35.0	26,300
23	787	424	-16.6	65,000	106	1769	509	-1.5	56,100	187	773	296	-17.0	90,800
24	313	417	<-35.0	66,000	107	1585	720	-3.6	42,500	188	1538	807	-4.2	38,400
25	807	516	-16.1	55,500	108	1692	807	-2.4	38,300	191	1560	674	-3.9	44,900
27	1184	524	-9.0	54,900	109	1482	593	-4.8	49,700	192	1818	687	-0.9	44,200
28	1263	446	-8.0	62,400	110	778	516	-16.9	55,500	193	1469	555	-5.0	52,400
29	743	605	-17.8	49,000	111	1728	700	-2.0	43,500	194	1380	266	-6.4	101,600
30	768	112	-17.2	348,600	113	1191	680	-8.9	44,500	195	784	632	-16.7	47,300
32	1216	417	-8.6	66,000	114	1298	185	-7.5	160,800	196	1227	1185	-8.4	23,700
33	1145	445	-9.5	62,500	115	682	907	-19.6	34,100	197	667	553	-20.1	52,600
34	1037	555	-11.3	52,400	116	1146	610	-9.5	48,700	198	2006	681	>0.0	44,500
35	863	412	-14.9	66,600	117	1548	849	-4.1	36,500	199	1711	674	-2.2	44,900
36	712	606	-18.7	48,900	118	1050	577	-11.1	50,800	200	872	424	-14.7	65,000
38	763	694	-17.3	43,800	120	1530	828	-4.3	37,400	201	292	435	<-35.0	63,700
39	304	470	<-35.0	59,800	121	638	423	-15.4	65,200	202	736	253	-18.0	107,800
41	1165	569	-9.2	51,400	122	1572	712	-3.8	42,900	203	786	829	-16.7	37,400
42	684	607	-19.6	48,800	123	23	1433	<-35.0	15,300	204	1224	589	-8.5	50,000
43	1318	589	-7.3	50,000	124	621	1474	-21.9	13,900	205	439	983	-30.9	31,100
44	1924	362	-0.1	74,600	125	1298	862	-7.5	36,000	206	1994	571	>0.0	51,300
46	1203	586	-8.7	50,200	126	872	921	-14.7	33,500	207	1895	687	-0.3	44,200
47	1391	447	-6.3	62,300	127	1000	717	-12.0	42,600	208	240	1418	<-35.0	15,800
48	309	454	<-35.0	61,500	128	1229	311	-8.4	86,100	210	1700	499	-2.3	57,000
49	605	587	-22.5	50,100	129	1422	832	-5.8	37,300	211	902	517	-14.1	55,400
50	621	535	-21.8	53,900	130	1776	499	-1.4	57,000	213	1087	684	-10.4	44,400
51	1113	522	-10.0	55,000	131	1930	757	-0.1	40,700	214	1340	668	-7.0	45,200
52	1820	499	-0.9	57,000	132	660	537	-20.4	53,800	215	1591	495	-3.5	57,300
53	725	177	-18.3	170,800	133	666	1019	-20.2	29,700	216	1585	755	-3.6	40,700
54	2001	500	>0.0	56,900	134	1271	862	-7.9	36,000	217	1159	393	-9.3	69,300
55	722	830	-18.4	37,300	135	1161	1389	-9.3	16,800	218	931	572	-13.5	51,200
56	678	533	-19.8	54,100	136	453	1063	-29.7	28,100	219	713	177	-18.7	170,500
57	1682	302	-2.5	89,000	137	1858	823	-0.6	37,700	220	1479	911	-4.9	33,900
58	1091	580	-10.3	50,600	138	1504	697	-4.6	43,700	221	965	927	-12.8	33,300
59	1171	585	-9.2	50,300	139	1488	707	-4.8	43,200	223	934	716	-13.5	42,700
60	1400	624	-6.2	47,800	140	1689	756	-2.4	40,700	225	1812	1045	-1.0	28,800
61	1853	508	-0.6	56,200	141	311	1417	<-35.0	15,800	226	821	411	-15.8	66,800
62	1888	567	-0.4	51,500	142	1366	915	-6.7	33,800	227	1586	1483	-3.6	13,600
65	735	297	-18.1	90,500	143	1429	346	-5.7	77,900	228	1065	567	-10.8	51,600
66	1263	312	-8.0	85,900	144	615	1017	-22.1	29,800	229	1577	890	-3.7	34,800
67	1252	407	-8.1	67,300	145	2006	566	>0.0	51,600	230	1458	496	-5.2	57,300
68	779	692	-16.8	43,900	146	2006	518	>0.0	55,300	232	1440	849	-5.5	36,500
69	1064	296	-10.8	90,800	147	1070	1108	-10.7	26,500	234	1692	489	-2.4	57,900
71	656	589	-20.6	50,000	148	1347	578	-6.9	50,800	235	618	1004	-22.0	30,300
72	638	545	-21.2	53,100	149	541	1481	-25.7	13,700	236	920	1138	-13.7	25,400
73	1582	583	-3.6	50,400	150	1645	760	-2.8	40,500	237	952	1008	-13.1	30,200
74	1570	556	-3.8	52,300	151	1269	236	-7.9	117,000	238	1611	541	-3.2	53,500
75	1264	621	-8.0	48,000	152	1507	911	-4.5	33,900	239	1489	720	-4.8	42,500
76	1338	564	-7.0	51,800	153	1722	448	-2.1	62,100	240	501	448	-27.7	62,100
77	1833	363	-0.8	74,400	154	932	503	-13.5	56,600	241	1820	569	-0.9	51,400
78	1767	565	-1.5	51,700	155	1031	294	-11.4	91,400	242	1357	658	-6.8	45,800
79	925	738	-13.6	41,600	156	1970	684	>0.0	44,400	243	711	1182	-18.7	23,800
80	534	698	-26.1	43,600	157	1258	183	-8.1	162,400	244	1855	621	-0.6	48,000
81	1811	363	-1.0	74,500	158	1275	417	-7.8	65,900	245	1189	474	-8.9	59,300
82	1412	681	-6.0	44,500	159	1663	820	-2.6	37,800	246	551	459	-25.1	61,000
83	1471	347	-5.0	77,500	160	1034	527	-11.4	54,600	247	1348	604	-6.9	49,100
84	1662	563	-2.7	51,800	161	1953	771	>0.0	40,000	248	460	448	-29.3	62,100
85	1596	479	-3.4	58,900	162	1020	1482	-11.6	13,700	249	1733	451	-1.9	61,800
86	1817	301	-0.9	89,100	164	1566	806	-3.8	38,400	250	1974	788	>0.0	39,200
87	516	1371	-27.0	17,400	166	1905	565	-0.2	51,700	251	808	392	-16.1	69,500
88	1589	698	-3.5	43,600	167	1340	181	-7.0	164,900	252	874	553	-14.6	52,500
89	1706	719	-2.2	42,500	168	1506	583	-4.6	50,400	253	753	848	-17.6	36,500
90	651	329	-20.8	81,700	169	1338	678	-7.0	44,700	254	995	450	-12.1	61,900
91	1415	710	-6.0	43,000	170	1969	541	>0.0	53,500	255	1690	679	-2.4	44,600
92	1773	545	-1.4	53,200	171	800	378	-16.3	71,800	256	994	1006	-12.1	30,200
93	1338	446	-7.0	62,300	172	476	958	-28.7	32,100	257	508	464	-27.4	60,400
94	1708	696	-2.2	43,700	173	919	1314	-13.7	19,300	258	1517	820	-4.4	37,800

<sup>a)</sup> Master table of proteins in the rat liver database, showing spot master number, gel position (x and y), isoelectric point relative to CPK standards, and predicted molecular mass (from the standard curve of Fig. 8).

MSN	X	Y	CPKd	SDSMW	MSN	X	Y	CPKd	SDSMW	MSN	X	Y	CPKd	SDSMW
259	1796	961	-1.1	31,900	345	1006	578	-11.9	50,800	426	1296	704	-7.8	43,300
260	661	1361	-20.4	17,700	346	1095	640	-10.3	46,800	427	810	843	-16.0	36,800
261	1725	679	-2.0	44,600	347	625	728	-21.7	42,000	428	1565	303	-3.9	88,700
262	496	1127	-28.0	25,800	348	361	983	-35.3	31,100	429	1259	847	-8.0	36,600
263	1063	172	-10.9	177,400	349	110	1343	<-35.0	18,300	430	1253	562	-8.1	51,900
265	1390	673	-6.3	45,000	350	521	1130	-26.7	25,700	431	734	1426	-18.1	15,500
266	510	437	-27.3	63,400	351	912	619	-13.9	48,100	432	483	433	-28.5	63,900
267	660	1038	-20.4	29,000	352	1574	530	-3.7	54,300	434	518	1041	-26.9	26,900
268	430	961	-31.0	31,900	353	961	912	-12.9	33,900	435	1020	1170	-11.6	24,300
269	1044	606	-11.2	48,900	354	706	762	-18.9	40,400	436	1122	196	-9.8	147,600
270	2019	853	>0.0	36,300	355	1450	830	-5.3	37,300	437	1870	673	-0.5	45,000
271	857	422	-15.0	65,200	356	1374	1152	-6.5	24,900	438	435	1102	-31.0	26,700
272	895	968	-14.2	31,700	357	474	967	-28.7	30,600	439	86	847	<-35.0	36,600
274	1292	712	-7.6	42,900	358	798	346	-16.3	77,800	440	1740	544	-1.8	53,200
275	1350	580	-6.9	49,900	359	764	338	-17.3	79,400	441	599	1571	-22.8	10,800
276	1670	1089	-2.6	27,100	360	1384	1068	-6.4	27,900	443	743	335	-17.8	80,100
277	688	538	-19.4	53,700	361	1713	769	-2.1	40,100	446	801	668	-16.2	45,200
278	961	718	-13.0	42,600	362	1161	859	-9.3	36,100	447	1050	926	-11.1	33,300
279	879	570	-14.5	51,300	363	914	1156	-13.8	24,800	448	1245	1298	-8.2	19,800
281	1848	1084	-0.7	27,300	364	412	435	-32.0	63,700	449	1576	1516	-3.7	12,600
282	1505	525	-4.6	54,800	365	741	486	-17.9	58,200	450	1818	1021	-0.9	29,600
283	1313	1147	-7.3	25,100	366	878	1503	-14.6	13,000	451	1094	440	-10.3	63,100
284	1314	829	-7.3	37,400	367	1560	935	-3.9	33,000	452	1945	802	>0.0	38,600
285	1332	408	-7.1	67,200	368	983	520	-12.4	55,200	453	1652	894	-2.8	34,600
286	1277	652	-7.8	46,100	369	434	441	-31.0	63,000	454	1403	500	-6.1	56,900
288	1391	824	-6.3	37,600	370	639	610	-21.2	48,700	456	1394	718	-6.3	42,600
289	1147	579	-9.5	50,700	371	1587	860	-3.6	36,100	457	905	436	-14.0	63,500
290	925	511	-13.6	55,900	372	1875	762	-0.5	40,400	459	1038	581	-11.3	50,500
291	787	1476	-16.6	13,900	373	1351	1059	-6.8	28,300	460	1598	294	-3.4	91,400
292	1462	818	-5.1	37,800	374	1506	715	-4.6	42,700	461	1528	863	-4.3	35,900
293	531	449	-26.3	62,000	375	1823	532	-0.9	54,200	462	1098	1137	-10.2	25,400
294	860	698	-14.9	43,600	376	254	417	<-35.0	65,900	463	849	1125	-15.2	25,800
295	1162	609	-9.3	48,700	377	1409	583	-6.1	50,400	464	1814	1072	-0.9	27,800
296	218	814	<-35.0	38,000	378	621	494	-21.8	57,500	465	1388	481	-6.3	58,700
297	1377	979	-6.5	31,300	379	1017	595	-11.7	49,600	466	1194	1084	-8.9	27,300
299	913	1523	-13.9	12,400	381	953	598	-13.1	49,400	468	577	467	-23.9	60,100
300	2012	667	>0.0	45,300	382	856	674	-15.0	44,900	469	1140	888	-9.6	34,900
301	702	178	-19.0	169,200	383	1252	258	-8.1	105,300	470	1797	524	-1.1	54,800
302	494	1280	-28.1	20,400	384	1699	1518	-2.3	12,500	471	1293	1133	-7.6	25,500
303	403	1008	-32.6	30,100	385	1042	493	-11.2	57,500	472	618	655	-21.9	46,000
304	1843	1585	-0.7	10,300	386	1490	583	-4.7	50,400	473	2009	299	>0.0	89,900
305	1049	583	-11.1	49,800	387	1554	603	-4.0	49,100	474	1205	215	-8.7	131,300
306	1608	989	-3.3	30,900	388	1193	404	-8.9	67,700	475	1035	788	-11.4	39,200
307	1219	916	-8.5	33,700	389	1374	902	-6.5	34,300	476	160	155	<-35.0	207,600
308	1627	755	-3.0	40,700	390	1456	969	-5.2	31,700	477	469	1370	-28.9	17,400
309	1524	892	-4.4	34,700	391	718	690	-18.5	44,000	478	599	662	-22.8	45,600
310	1769	1028	-1.5	29,400	392	1799	732	-1.1	41,900	479	1009	540	-11.8	53,500
311	1609	1451	-3.3	14,700	393	1482	758	-4.8	40,600	480	1216	235	-8.6	117,400
312	266	1408	<-35.0	16,100	394	1227	1461	-8.4	14,400	482	816	346	-15.9	77,800
313	1902	1365	-0.3	17,600	395	1530	577	-4.3	50,800	483	693	673	-19.3	44,900
314	1316	1395	-7.3	16,600	396	1410	755	-6.0	40,800	485	1608	1013	-3.3	30,000
315	1341	523	-7.0	54,900	397	912	256	-13.9	106,400	486	478	599	-28.6	49,300
318	1104	1053	-10.1	28,500	399	1465	1063	-5.0	28,100	487	1025	607	-11.5	48,800
320	1480	1459	-4.9	14,400	400	1473	450	-4.9	61,900	488	1045	1186	-11.2	23,700
321	850	603	-15.1	49,100	401	1029	1140	-11.5	25,300	489	1609	301	-3.3	89,200
322	1454	1494	-5.3	13,300	403	1516	754	-4.4	40,800	490	775	1289	-17.0	20,100
323	670	626	-20.0	47,700	404	1495	554	-4.7	52,500	491	692	178	-19.3	169,300
324	655	101	-20.6	420,500	405	1525	1092	-4.3	27,100	492	1100	964	-10.2	31,800
325	1521	675	-4.4	44,800	406	723	252	-18.4	108,000	493	1760	776	-1.6	39,700
326	1587	677	-3.6	44,700	409	650	663	-20.8	45,500	494	882	247	-14.5	110,700
327	1388	409	-6.3	67,000	410	1501	478	-4.6	59,000	495	470	1258	-28.9	21,200
328	448	1291	-30.0	20,100	411	936	1057	-13.4	28,300	496	494	1436	-28.1	15,200
330	1608	751	-3.3	40,900	412	350	1120	-35.9	26,000	497	980	852	-12.5	36,400
331	1566	697	-3.8	43,700	413	1033	538	-11.4	53,700	499	1414	546	-6.0	53,100
332	531	471	-26.3	59,600	415	737	425	-18.0	64,900	500	1234	1072	-8.3	27,800
333	784	1156	-16.7	24,700	416	1578	606	-3.7	48,900	501	1246	659	-8.2	45,700
334	1059	407	-10.9	67,300	417	646	496	-21.0	57,300	502	824	792	-15.7	39,000
335	1593	303	-3.5	88,500	418	1695	482	-2.3	58,600	503	1246	1134	-8.2	25,500
336	1616	598	-3.2	49,400	419	725	770	-18.3	40,000	504	1115	1407	-9.9	16,200
338	1854	1004	-0.6	30,300	420	1289	1041	-7.7	28,900	505	1189	391	-8.9	69,700
339	1265	888	-8.0	34,900	421	1171	912	-9.1	33,900	506	1578	402	-3.7	68,000
340	581	585	-23.6	50,300	422	599	162	-22.8	193,700	507	787	250	-16.6	109,000
341	1497	1047	-4.7	28,700	423	929	856	-13.6	36,200	508	979	552	-12.5	52,600
343	1351	265	-6.8	102,200	424	739	625	-17.9	47,700	509	1153	619	-9.4	48,100
344	1813	549	-0.9	52,800	425	1490	965	-4.7	31,800	510	1730	1006	-2.0	30,200

MSN	X	Y	CPKbl	SDSMW	MSN	X	Y	CPKbl	SDSMW	MSN	X	Y	CPKbl	SDSMW
511	809	484	-16.0	58,400	596	619	269	-21.9	100,500	674	1661	448	-2.7	62,100
512	1099	533	-10.2	54,100	597	1176	461	-9.1	60,700	675	1523	562	-4.4	51,900
513	1696	1034	-2.3	29,200	598	1465	1044	-5.0	28,800	676	708	642	-18.8	46,700
514	948	636	-13.2	47,100	599	741	1188	-17.9	23,600	677	919	615	-13.7	48,300
515	481	543	-28.5	53,400	600	907	402	-14.0	68,000	678	1085	551	-10.5	52,700
516	1334	1044	-7.1	28,800	601	687	658	-19.5	45,800	679	600	923	-22.7	33,400
517	868	1021	-14.8	29,700	602	712	1138	-18.7	25,400	680	1237	1004	-8.3	30,300
518	798	779	-16.3	39,600	603	898	181	-14.1	165,200	681	1103	283	-10.1	95,100
519	822	670	-15.7	45,100	604	783	1461	-16.7	14,400	682	1406	477	-6.1	59,100
520	632	165	-21.5	189,000	605	736	223	-18.0	125,300	683	1596	249	-3.4	109,800
521	1332	830	-7.1	37,300	606	629	273	-21.6	98,700	684	555	699	-24.8	43,500
522	603	1104	-22.6	26,600	607	1064	286	-10.8	94,000	685	1167	1313	-9.2	19,300
523	1190	309	-8.9	86,800	608	883	503	-14.5	56,700	686	1932	790	0.0	39,100
524	479	1226	-28.6	22,300	609	2012	610	>0.0	48,700	687	1545	619	-4.1	48,100
525	768	1066	-17.2	28,000	610	1255	903	-8.1	34,200	688	1456	764	-5.2	40,300
526	747	1016	-17.7	29,800	612	1103	391	-10.1	69,600	689	1011	953	-11.8	32,300
527	1170	231	-9.2	119,600	613	778	265	-16.9	102,000	690	1995	270	>0.0	100,200
528	1502	542	-4.6	53,400	614	824	518	-15.7	55,400	691	812	888	-16.0	34,900
530	1728	620	-2.0	48,000	615	1095	195	-10.3	149,100	692	1154	1461	-9.4	14,400
532	507	1011	-27.4	30,000	616	1759	478	-1.6	59,000	693	1993	819	>0.0	37,800
533	870	489	-14.7	57,900	617	994	372	-12.1	72,900	694	1628	656	-3.0	45,900
534	1347	1085	-6.9	27,300	618	751	374	-17.6	72,400	695	928	254	-13.6	107,000
535	1513	346	-4.5	77,800	619	1429	518	-5.7	55,300	696	1854	715	-0.6	42,700
536	308	654	<-35.0	46,000	620	1050	520	-11.1	55,200	697	1997	345	>0.0	78,000
538	1851	689	-0.7	44,100	621	923	1105	-13.7	26,600	698	957	563	-13.0	51,800
539	1463	982	-5.1	31,100	622	1462	622	-5.1	47,900	699	1540	730	-4.2	42,000
540	909	561	-13.9	52,000	623	759	225	-17.4	124,000	702	577	900	-23.8	34,400
541	625	289	-21.7	93,100	624	758	1038	-17.4	29,000	703	1610	562	-3.2	51,900
542	1164	198	-9.2	146,200	625	1438	606	-5.5	48,900	705	1278	571	-7.8	51,200
543	803	655	-16.2	45,900	626	1096	1089	-10.2	27,200	706	1841	704	-0.7	43,300
544	1259	1143	-8.0	25,200	627	942	548	-13.3	53,000	707	1018	1386	-11.7	16,900
545	856	1526	-15.0	12,200	628	809	621	-16.0	48,000	709	1074	1145	-10.7	25,100
546	803	1071	-16.2	27,800	629	899	979	-14.1	31,300	710	293	889	<-35.0	34,800
547	1162	274	-9.3	98,400	630	1135	1321	-9.6	19,100	712	720	412	-18.5	66,600
548	128	1321	<-35.0	19,000	631	979	615	-12.5	48,300	713	1386	841	-6.4	36,800
549	1355	1122	-6.8	25,900	632	1542	1076	-4.1	27,600	714	1328	263	-7.1	103,100
550	595	866	-23.0	35,800	633	1345	814	-6.9	38,000	715	698	433	-19.1	63,900
552	1369	494	-6.6	57,500	634	409	950	-32.2	32,400	716	701	481	-19.0	58,700
553	992	405	-12.2	67,600	635	1165	704	-9.2	43,300	717	1875	699	-0.5	43,600
555	1125	410	-9.8	66,900	636	774	604	-17.0	49,000	718	575	702	-23.9	43,400
556	705	975	-18.9	31,400	637	1263	524	-8.0	54,800	719	1216	204	-8.6	140,400
557	1477	1030	-4.9	29,300	638	952	411	-13.1	66,700	721	1069	464	-10.8	60,400
558	980	583	-12.5	50,400	639	1717	575	-2.1	51,000	722	1272	506	-7.9	56,400
559	700	1109	-19.1	26,400	640	994	292	-12.1	92,000	723	958	822	-13.0	37,700
560	1028	621	-11.5	48,000	641	165	1224	<-35.0	22,400	724	763	395	-17.3	69,100
562	898	794	-14.1	38,900	642	803	251	-16.2	108,900	725	720	916	-18.5	33,700
564	789	1446	-16.6	14,900	643	719	296	-18.5	90,700	726	1476	415	-4.9	66,200
565	777	766	-16.9	40,200	644	1100	294	-10.2	91,400	727	1846	473	-0.7	59,400
566	980	328	-12.5	81,900	645	534	1263	-26.1	21,000	728	510	783	-27.3	39,400
567	1519	611	-4.4	48,600	646	1153	1038	-9.4	29,000	729	1217	1126	-8.6	25,800
569	1212	661	-8.6	45,600	648	1246	204	-8.2	140,000	730	1858	724	-0.6	42,300
570	760	594	-17.4	49,700	649	14	1406	<-35.0	16,200	731	665	765	-20.2	40,300
571	618	956	-21.9	32,100	650	1713	1049	-2.1	28,600	733	1321	312	-7.2	85,900
573	1142	771	-9.6	40,000	651	1986	1183	>0.0	23,800	734	719	427	-18.5	64,600
574	532	787	-26.2	39,300	652	1378	816	-6.5	38,000	735	1101	473	-10.2	59,500
575	771	250	-17.1	109,200	653	1442	1165	-5.5	24,400	736	1359	569	-6.7	51,400
576	1068	534	-10.8	54,100	654	650	806	-20.8	38,400	738	696	220	-19.2	127,600
577	822	734	-15.7	41,800	655	1111	551	-10.0	52,700	739	687	409	-19.5	67,000
578	914	754	-13.8	40,800	656	1095	861	-10.3	36,000	740	1205	256	-8.7	106,200
579	1064	794	-10.8	38,900	657	1524	540	-4.4	53,600	741	995	563	-12.1	51,900
580	1524	714	-4.4	42,800	658	1777	860	-1.4	36,000	742	898	596	-14.1	49,500
581	1392	783	-6.3	39,400	659	391	584	-33.4	50,400	743	881	181	-14.5	165,900
582	982	686	-12.4	44,200	660	977	565	-12.5	51,700	744	1951	686	>0.0	44,200
584	1487	672	-4.8	45,000	661	658	166	-20.5	187,500	745	726	168	-18.3	183,600
585	758	731	-17.4	41,900	662	732	312	-18.1	86,100	746	999	643	-12.0	46,600
586	687	1152	-19.5	24,900	663	1787	567	-1.2	51,500	748	182	1503	<-35.0	13,000
587	930	523	-13.5	55,000	664	888	268	-14.4	100,900	749	2005	649	>0.0	46,300
588	1888	774	-0.4	39,900	665	889	775	-14.3	39,800	750	1448	575	-5.4	51,000
589	642	485	-21.1	58,300	666	715	221	-18.6	126,300	751	792	266	-16.5	101,900
590	1317	519	-7.3	55,300	667	781	227	-16.8	122,400	752	469	296	-28.9	90,600
591	65	1548	<-35.0	11,500	668	646	165	-21.0	189,100	754	664	254	-20.3	107,000
592	1014	614	-11.7	48,400	669	1116	353	-9.9	76,300	755	1195	184	-8.8	161,000
593	732	176	-18.1	172,300	670	1382	643	-6.4	46,600	756	1821	1113	-0.9	26,300
594	1627	478	-3.0	59,000	671	547	789	-25.3	39,200	757	909	246	-13.9	111,000
595	1009	1426	-11.8	15,500	673	984	746	-12.4	41,200	760	790	133	-16.5	264,900

MSN	X	Y	CPKdI	SDSMW	MSN	X	Y	CPKdI	SDSMW	MSN	X	Y	CPKdI	SDSMW
761	1399	733	-6.2	41,800	648	1863	271	-0.6	99,500	939	1197	827	-8.8	37,500
763	1416	1085	-5.9	27,300	649	1166	523	-9.2	54,900	941	1765	885	-1.5	35,000
764	2020	569	>0.0	51,400	650	1535	1024	-4.2	29,600	942	602	472	-22.7	59,600
765	651	475	-20.8	59,300	851	1035	826	-11.4	37,500	943	312	488	<-35.0	57,100
766	1052	1149	-11.1	25,000	852	834	542	-15.5	53,400	944	993	491	-12.1	57,700
767	1968	468	>0.0	59,900	855	499	220	-27.8	127,100	945	1300	269	-7.5	100,300
768	1330	685	-7.1	44,300	856	1063	194	-10.9	150,500	946	630	423	-21.6	65,100
769	1970	613	>0.0	48,500	857	887	890	-14.4	34,800	947	187	736	<-35.0	41,600
770	857	617	-15.0	48,200	858	1448	639	-5.4	46,900	948	1380	344	-6.5	78,200
771	1337	974	-7.0	31,500	859	706	311	-18.9	86,200	949	1766	665	-1.5	45,400
773	1576	502	-3.7	56,700	860	1070	1066	-10.7	28,000	950	1038	193	-11.3	151,000
775	969	824	-12.8	37,600	861	472	347	-28.8	77,600	951	860	152	-14.9	213,000
776	1438	708	-5.5	43,100	862	674	480	-19.9	58,800	952	957	701	-13.0	43,400
777	1539	458	-4.2	61,000	864	1307	499	-7.4	57,000	954	503	547	-27.6	53,000
778	850	434	-15.1	63,800	865	645	887	-21.0	34,900	955	1938	712	>0.0	42,900
779	700	411	-19.1	66,800	866	827	1004	-15.6	30,300	957	1010	816	-11.8	37,900
780	1052	1136	-11.1	25,500	868	685	494	-19.5	57,400	959	768	174	-17.2	174,900
784	1413	529	-6.0	54,400	869	1807	402	-1.0	68,000	960	596	419	-23.0	65,700
785	1364	885	-6.7	35,000	870	1323	783	-7.2	39,400	961	557	409	-24.8	67,100
786	1822	835	-0.9	37,100	871	1228	1031	-8.4	29,300	962	887	320	-14.4	83,900
787	893	392	-14.3	68,500	872	1904	346	-0.3	77,700	963	564	334	-24.5	80,500
790	616	882	-22.0	35,100	873	556	647	-24.8	46,400	964	969	1155	-12.8	24,800
791	451	1429	-29.8	15,400	874	1540	756	-4.2	40,700	965	671	255	-20.0	106,600
792	777	377	-16.9	72,000	875	1566	777	-3.8	39,700	966	1204	798	-8.7	38,700
793	1536	1543	-4.2	11,700	876	1198	351	-8.8	76,800	967	910	154	-13.9	210,300
794	1461	807	-5.1	38,300	877	1076	720	-10.6	42,500	968	609	1048	-22.3	28,700
796	388	546	-33.6	53,100	878	1161	1111	-9.3	26,400	969	1285	206	-7.7	138,900
797	1126	212	-9.8	133,700	879	647	757	-20.9	40,700	970	822	232	-15.8	119,300
798	933	437	-13.5	63,400	880	1756	594	-1.6	49,700	971	976	437	-12.6	63,400
799	1420	593	-5.9	49,800	881	1543	278	-4.1	97,100	972	403	567	-32.6	51,600
800	1759	279	-1.6	96,500	883	1432	890	-5.7	34,800	974	279	495	<-35.0	57,400
801	624	865	-21.7	35,800	884	922	689	-13.7	44,100	975	844	981	-15.3	31,200
802	898	547	-14.2	53,000	885	1103	414	-10.1	66,400	976	1124	295	-9.8	91,100
803	1775	1468	-1.4	14,200	886	1501	607	-4.6	48,900	977	994	664	-12.1	45,400
804	573	196	-24.0	148,400	887	798	1103	-16.3	26,600	978	1612	642	-3.2	46,700
805	203	494	<-35.0	57,400	888	636	634	-21.3	47,200	979	749	1141	-17.7	25,300
806	980	1039	-12.5	29,000	889	951	759	-13.1	40,600	980	1064	642	-10.8	46,700
807	902	308	-14.1	87,200	890	717	548	-18.6	52,900	981	1197	911	-8.8	33,900
808	625	827	-21.7	37,500	891	1123	229	-9.8	121,200	983	1762	1508	-1.6	12,800
809	1851	1015	-0.7	29,900	892	891	413	-14.3	66,400	984	1344	317	-6.9	84,700
810	440	573	-30.9	51,100	894	1245	234	-8.2	117,800	985	1024	1105	-11.5	26,600
811	1358	249	-6.8	109,700	895	1962	346	>0.0	77,700	987	739	1159	-17.9	24,600
812	851	393	-15.1	69,400	896	1322	626	-7.2	47,700	988	816	555	-15.9	52,400
813	745	1246	-17.8	21,600	897	420	570	-31.4	51,300	990	785	361	-16.7	74,900
814	2028	810	>0.0	38,200	898	662	428	-20.3	64,500	991	1159	317	-9.3	84,500
815	1086	645	-10.4	46,500	899	845	243	-15.3	113,000	992	1090	928	-10.4	33,300
816	629	313	-21.6	85,700	900	624	703	-21.7	43,400	993	1030	701	-11.5	43,400
817	1376	1177	-6.5	24,000	901	931	1094	-13.5	27,000	994	847	811	-15.2	38,200
818	1771	790	-1.4	39,100	903	799	229	-16.3	121,000	995	902	461	-14.1	60,700
819	1045	263	-11.2	103,100	904	765	520	-17.2	55,200	996	888	847	-14.4	36,600
820	984	362	-12.4	74,600	905	775	889	-17.0	34,800	997	1815	579	-0.9	50,700
821	1712	279	-2.2	96,700	907	888	824	-14.4	37,600	998	1205	504	-8.7	56,500
822	1256	205	-8.1	139,200	908	828	1303	-15.6	19,700	999	617	289	-22.0	93,100
823	1517	654	-4.4	46,000	910	681	1544	-19.7	11,700	1000	968	290	-12.8	92,700
824	1442	449	-5.5	62,000	911	1544	301	-4.1	89,100	1001	970	771	-12.7	40,000
825	1240	513	-8.3	55,800	913	1606	387	-3.3	70,400	1002	1736	478	-1.9	58,900
826	1309	1014	-7.4	29,900	914	1237	688	-8.3	44,100	1003	643	1184	-21.1	23,700
827	2012	708	>0.0	43,100	916	1442	749	-5.5	41,100	1006	822	487	-15.8	58,100
828	937	1405	-13.4	16,200	917	1260	367	-8.0	73,700	1007	875	279	-14.6	96,400
830	1342	756	-7.0	40,700	919	764	1541	-17.3	11,700	1009	291	644	<-35.0	46,600
831	562	826	-24.5	37,500	920	1133	1123	-9.7	25,900	1010	1386	745	-6.4	41,200
832	1073	1039	-10.7	29,000	921	1123	380	-9.8	71,500	1011	459	541	-29.4	53,500
833	481	820	-28.5	37,800	923	829	242	-15.6	113,200	1012	679	661	-19.7	45,600
834	501	581	-27.8	50,500	924	1131	318	-9.7	84,300	1013	1818	1128	-0.9	25,800
837	751	748	-17.6	41,100	925	1441	874	-5.5	35,400	1014	1032	634	-11.4	47,200
838	635	833	-21.3	37,200	926	679	219	-19.7	128,200	1015	1629	994	-3.0	30,700
839	1494	459	-4.7	60,900	927	1487	1191	-4.8	23,500	1016	1311	1134	-7.4	25,500
840	1952	301	>0.0	89,300	928	1082	775	-10.5	39,800	1017	1722	424	-2.0	65,000
841	1585	1080	-3.6	27,500	929	1231	816	-8.4	38,000	1018	1015	743	-11.7	41,300
842	571	1312	-24.1	19,400	931	1609	670	-3.3	45,100	1020	1574	1219	-3.7	22,500
843	1325	649	-7.2	46,300	932	810	900	-16.0	34,400	1021	781	484	-16.8	58,400
844	1727	301	-2.0	89,200	933	965	520	-12.8	55,100	1022	1129	83	-9.7	591,300
845	630	679	-21.5	44,600	934	947	462	-13.2	60,600	1023	812	317	-15.9	84,600
846	2016	905	>0.0	34,200	936	865	843	-14.8	36,800	1024	785	446	-16.7	62,400
847	673	1200	-19.9	23,200	937	1421	1056	-5.9	28,400	1025	1290	739	-7.7	41,500

MSN	X	Y	CPKpl	SDSMW	MSN	X	Y	CPKpl	SDSMW	MSN	X	Y	CPKpl	SDSMW
1026	405	552	-32.5	52,600	1153	921	1158	-13.7	24,700	1246	547	577	-25.3	50,800
1027	1298	848	-7.5	36,500	1154	1564	864	-3.5	35,900	1247	530	576	-26.3	50,900
1028	856	547	-15.0	53,000	1161	637	400	-21.3	68,400	1249	516	572	-27.0	51,200
1030	1284	226	-7.7	123,200	1162	623	397	-21.8	68,800	1250	973	536	-12.7	53,900
1031	886	822	-12.3	37,700	1163	665	397	-20.2	68,700	1251	607	532	-22.4	54,200
1032	1547	403	-4.1	67,900	1168	564	528	-24.4	54,500	1252	665	529	-20.2	54,400
1033	1381	551	-6.4	52,700	1170	552	529	-25.0	54,500	1253	899	766	-14.1	40,200
1034	1525	486	-4.3	57,200	1171	538	524	-25.9	54,800	1254	1311	746	-7.4	41,200
1035	1128	645	-9.7	46,500	1172	545	514	-25.5	55,700	1255	1300	761	-7.5	40,400
1036	1226	274	-8.5	98,300	1174	1099	522	-10.2	55,000	1257	1938	712	0.0	42,900
1039	1761	262	-1.6	103,600	1176	1304	586	-7.5	50,200	1258	1806	718	-1.0	42,600
1040	541	839	-25.7	36,900	1177	1366	539	-6.6	53,700	1259	1727	715	-2.0	42,700
1041	818	910	-15.8	34,000	1178	1608	702	-3.3	43,400	1260	1629	713	-3.0	42,800
1044	1036	485	-11.3	58,300	1179	1485	224	-4.8	124,900	1261	1555	717	-4.0	42,600
1045	1439	407	-5.5	67,300	1180	1459	224	-5.2	124,900	1262	1468	717	-5.0	42,600
1047	1540	250	-4.2	109,200	1181	1431	223	-5.7	125,100	1263	1413	722	-6.0	42,400
1048	1576	635	-3.7	47,100	1182	1407	223	-6.1	125,200	1264	1340	717	-7.0	42,600
1049	1089	411	-10.4	66,700	1183	1383	224	-6.4	124,700	1265	1263	717	-8.0	42,600
1050	849	1040	-13.2	28,900	1184	1454	182	-5.3	164,400	1266	1182	720	-9.0	42,500
1051	426	818	-31.1	37,800	1185	1422	183	-5.8	162,600	1267	1110	717	-10.0	42,600
1052	1583	1385	-3.6	16,900	1186	1394	182	-6.3	164,300	1268	1055	717	-11.0	42,600
1053	779	1092	-16.8	27,000	1189	1171	214	-9.2	131,800	1269	999	717	-12.0	42,600
1054	1613	620	-3.2	48,000	1190	1457	286	-5.2	94,200	1270	959	715	-13.0	42,700
1055	1380	377	-6.5	72,000	1191	686	1114	-19.5	26,200	1271	905	712	-14.0	42,900
1056	284	663	<-35.0	45,500	1192	265	893	<-35.0	34,700	1272	857	714	-15.0	42,800
1058	1261	746	-8.0	41,200	1193	403	1292	-32.6	20,000	1273	810	705	-16.0	43,300
1060	393	605	-33.3	49,000	1194	344	1275	<-35.0	20,600	1274	774	711	-17.0	42,900
1061	1817	645	-0.9	46,600	1195	505	1311	-27.6	19,400	1277	737	708	-18.0	43,100
1062	1245	746	-8.2	41,200	1196	572	1293	-24.1	20,000	1278	702	711	-19.0	42,900
1064	1258	792	-8.1	39,000	1197	639	1502	-21.2	13,000	1279	671	710	-20.0	43,000
1065	705	934	-18.9	33,000	1198	637	1402	-21.3	16,300	1280	645	710	-21.0	43,000
1066	1181	734	-9.0	41,800	1199	614	1407	-22.1	16,200	1281	617	707	-22.0	43,100
1067	529	658	-26.3	45,800	1200	637	1431	-21.3	15,400	1282	595	704	-23.0	43,300
1068	508	696	-27.4	43,700	1201	1095	1394	-10.3	16,600	1283	573	700	-24.0	43,500
1069	1898	604	-0.3	49,100	1202	1719	1545	-2.1	11,600	1284	552	695	-25.0	43,700
1071	873	609	-14.7	48,700	1203	791	668	-16.5	45,200	1285	536	694	-26.0	43,800
1073	1768	1128	-1.5	25,800	1204	964	1021	-12.9	29,700	1286	515	687	-27.0	44,200
1075	836	773	-15.4	39,900	1205	313	195	<-35.0	148,700	1287	496	683	-28.0	44,400
1076	1863	861	-0.6	36,000	1208	306	194	<-35.0	149,800	1288	467	669	-29.0	45,200
1078	826	566	-15.7	51,600	1209	320	197	<-35.0	147,400	1289	447	667	-30.9	45,300
1081	971	483	-12.7	58,500	1210	326	197	<-35.0	146,600	1290	427	655	-31.0	45,900
1083	1697	202	-2.3	142,300	1211	394	294	-33.2	91,400	1291	412	655	-32.0	45,900
1085	1157	794	-9.4	38,900	1212	402	294	-32.7	91,200	1292	397	652	-33.0	46,100
1090	620	910	-21.9	34,000	1214	386	294	-33.7	91,400	1293	381	654	-34.0	46,000
1092	1867	597	-0.5	49,500	1215	641	329	-21.2	81,600	1294	365	653	-35.0	46,100
1093	2019	894	>0.0	34,600	1216	660	329	-20.4	81,600	1295	348	653	<-35.0	46,100
1094	1546	538	-4.1	53,700	1217	914	266	-13.8	101,800					
1095	1545	477	-4.1	59,100	1218	873	245	-14.7	112,000					
1098	61	935	<-35.0	33,000	1219	970	372	-12.7	72,900					
1099	1954	237	>0.0	116,000	1220	1021	298	-11.6	90,100					
1101	588	1048	-23.3	28,600	1221	1392	205	-6.3	139,500					
1102	1050	667	-11.1	45,200	1222	1354	203	-6.8	141,800					
1103	457	797	-29.5	38,800	1223	1362	205	-6.7	139,500					
1105	1884	532	-0.4	54,200	1224	673	540	-19.9	53,600					
1106	1714	649	-2.1	46,300	1225	614	542	-22.1	53,400					
1107	1717	546	-2.1	53,100	1226	603	539	-22.6	53,600					
1108	1976	722	>0.0	42,400	1227	696	623	-19.2	47,800					
1111	547	1066	-25.3	28,000	1228	707	628	-18.9	47,500					
1112	1348	621	-6.9	48,000	1229	475	447	-28.7	62,300					
1115	1385	762	-6.4	40,400	1230	466	1282	-29.0	20,400					
1116	1078	816	-10.6	38,000	1231	759	1461	-17.4	14,400					
1117	975	787	-12.6	39,300	1232	1324	1170	-7.2	24,200					
1118	1202	933	-8.7	33,100	1233	1583	1005	-3.6	30,300					
1119	1022	1076	-11.6	27,600	1234	1865	809	-0.6	38,200					
1120	1905	616	-0.3	48,300	1235	1812	817	-1.0	37,900					
1121	1512	1301	-4.5	19,700	1236	1411	703	-6.0	43,400					
1122	1114	677	-9.9	44,700	1237	1392	682	-6.3	44,500					
1123	1464	452	-5.1	61,700	1238	794	410	-16.4	66,900					
1125	1048	857	-11.1	36,200	1239	769	407	-17.1	67,300					
1126	1122	802	-9.8	38,600	1240	740	406	-17.9	67,500					
1128	1722	892	-2.1	34,700	1241	743	511	-17.8	55,900					
1133	1098	825	-10.2	37,500	1242	713	510	-18.7	56,000					
1139	1830	569	-0.8	51,400	1243	682	509	-19.6	56,100					
1147	764	1182	-17.3	23,800	1244	663	504	-20.3	56,500					
1148	1968	724	>0.0	42,300	1245	565	582	-24.4	50,500					

Table 2. Table of some identified proteins

POP name	Protein name	MSN's	Basis for identification
IDS:3_ALPHA_HDDH	3- $\alpha$ -hydroxysteroid-dihydrodiol-dehydrogenase, an enzyme of steroid metabolism	137, 159	Pure protein and antibody provided by Dr. T.M. Panning, Department of Pharmacology, School of Medicine, University of Pennsylvania.
IDS:ACTIN_BETA	$\beta$ cellular actin, a cytoskeletal protein	38	Homologous position with respect to other mammalian systems
IDS:ACTIN_GAMMA	$\gamma$ cellular actin, a cytoskeletal protein	68	Homologous position with respect to other mammalian systems
IDS:ALBUMIN	Serum albumin, mature form.	21, 28, 33	Predominance in rat plasma
IDS:APO_A-I	Apo A-I plasma lipoprotein, mature form (tentative)	236, 463	Presence in rat plasma, regulation by some lipid-lowering drugs
IDS:CALMODULIN	Calmodulin, an acidic cytosolic calcium-binding protein	123, 649	Homologous position with respect to other mammalian systems
IDS:CATALASE	Catalase (peroxisomal)	54, 61, 106	Presence in purified peroxisomes, similarity in position to mouse catalase
IDS:CPKSPOTS	Spots contributed by the CPK charge standards (not rat liver proteins)	1257 - 1295	
IDS:CPS	Carbamoyl phosphate synthase	114, 157, 167, 174, 1184, 1185, 1186, 1222	Pure protein provided by Dr. Margaret Marshall, Department of Pharmacology, Medical School, University of Wisconsin - Madison.
IDS:CYTOCHROME_B5	Cytochrome b5	87, 477	Pure protein provided by Dr. Andrew Parkinson, Department of Pharmacology, Toxicology and Therapeutics, University of Kansas Medical Center
IDS:FABP-L	Liver fatty-acid binding protein	227	Pure protein provided by Dr. Nathan Bass, Department of Medicine, University of California School of Medicine, San Francisco
IDS:HMG-COA_SYNTHASE	Cytosolic HMG-CoA Synthase	133, 144, 235, 413	Antibody provided by Dr. Michael Greenspan, Merck Sharp & Dohme Research Laboratories, Rahway, NJ
IDS:LAMIN_B	Lamin B, a nuclear protein	415, 734	Homologous position with respect to other mammalian systems
IDS:MITCON:1	Mitcon:1 (F1 ATPase $\beta$ subunit), a mitochondrial inner membrane protein equivalent to E.	17, 49, 71, 340, 1245, 1246, 1247, 1249	Homologous position with respect to other mammalian systems, presence in mitochondria
IDS:MITCON:2	Mitcon:2, a mitochondrial matrix stress protein	15, 25, 110, 1241, 1242, 1243, 1244	Homologous position with respect to other mammalian systems, presence in mitochondria
IDS:MITCON:3	Mitcon:3, a mitochondrial matrix stress protein, likely analog of NADPH cytochrome P-450 reductase, frequently co-induced with P-450's	18, 35, 226, 600, 1238, 1239, 1240	Homologous position with respect to other mammalian systems, presence in mitochondria
IDS:NADPH_P450_RED		175, 251, 812	Pure protein provided by Dr. Andrew Parkinson, Department of Pharmacology, Toxicology and Therapeutics, University of Kansas Medical Center
IDS:PDI	Protein disulphide isomerase 1	168, 1170, 1171, 1172	Sequence information obtained by R.M. Van Frank, Lilly Research Laboratories, Indianapolis
IDS:PLASMA_PROTEINS	Rat plasma proteins observed in liver	21, 28, 33, 44, 72, 102, 115, 197, 236, 246, 248, 257, 293, 332, 347, 364, 369, 419, 432, 463, 468, 518, 562, 605, 623, 666, 667, 725, 738, 790, 865, 903, 926	Plasma coelectrophoresis studies
IDS:PRO-ALBUMIN	Serum albumin precursor	47, 93	Relative position to mature albumin, presence in microsomes
IDS:PYRCARBOX	Pyruvate carboxylase	179, 1180, 1181, 1182, 1183	Pavlica, R.J., et al., BBA (1990) 1022 115-125.
IDS:SOD	Superoxide dismutase	135	Sequence information obtained by R.M. Van Frank, Lilly Research Laboratories, Indianapolis
IDS:TUBULIN_ALPHA	$\alpha$ tubulin, a cytoskeletal protein	56, 132, 1224, 1252	Homologous position with respect to other mammalian systems
IDS:TUBULIN_BETA	$\beta$ tubulin, a cytoskeletal protein	50, 1225, 1226, 1251	Homologous position with respect to other mammalian systems

Table 3. Computed pI's of two sets of carbamylated protein standards: Rabbit muscle CPK and human hemoglobin (Hb)

	Protein Name	PIR Name	#ASP 3.9	#GLU 4.1	#HIS 6.0	#LYS 10.8	#ARG 12.5	NH2- 7.0	Calc pI	Real CPK
0	Rabbit muscle CPK	KIRBCM	28	27	17	34	18	1	6.84	0.0
-1			28	27	17	33	18	1	6.67	-1
-2			28	27	17	32	18	1	6.54	-2
-3			28	27	17	31	18	1	6.42	-3
-4			28	27	17	30	18	1	6.31	-4
-5			28	27	17	29	18	1	6.21	-5
-6			28	27	17	28	18	1	6.12	-6
-7			28	27	17	27	18	1	6.03	-7
-8			28	27	17	26	18	1	5.94	-8
-9			28	27	17	25	18	1	5.85	-9
-10			28	27	17	24	18	1	5.76	-10
-11			28	27	17	23	18	1	5.67	-11
-12			28	27	17	22	18	1	5.58	-12
-13			28	27	17	21	18	1	5.48	-13
-14			28	27	17	20	18	1	5.39	-14
-15			28	27	17	19	18	1	5.29	-15
-16			28	27	17	18	18	1	5.20	-16
-17			28	27	17	17	18	1	5.12	-17
-18			28	27	17	16	18	1	5.04	-18
-19			28	27	17	15	18	1	4.96	-19
-20			28	27	17	14	18	1	4.89	-20
-21			28	27	17	13	18	1	4.83	-21
-22			28	27	17	12	18	1	4.77	-22
-23			28	27	17	11	18	1	4.71	-23
-24			28	27	17	10	18	1	4.66	-24
-25			28	27	17	9	18	1	4.61	-25
-26			28	27	17	8	18	1	4.56	-26
-27			28	27	17	7	18	1	4.52	-27
-28			28	27	17	6	18	1	4.48	-28
-29			28	27	17	5	18	1	4.44	-29
-30			28	27	17	4	18	1	4.40	-30
-31			28	27	17	3	18	1	4.36	-31
-32			28	27	17	2	18	1	4.32	-32
-33			28	27	17	1	18	1	4.29	-33
-34			28	27	17	0	18	1	4.25	-34
-35			28	27	17	0	18	0	4.22	-35
0	Hb-beta, human	HBHU	7	8	9	11	3	1	7.18	
-1			7	8	9	10	3	1	6.79	
-2			7	8	9	9	3	1	6.53	-1.8
-3			7	8	9	8	3	1	6.32	-3.2
-4			7	8	9	7	3	1	6.13	-5.3
-5			7	8	9	6	3	1	5.96	-7.2
-6			7	8	9	5	3	1	5.78	-10.0
-7			7	8	9	4	3	1	5.59	-12.3
-8			7	8	9	3	3	1	5.37	-15.5
-9			7	8	9	2	3	1	5.14	-18.0
-10			7	8	9	1	3	1	4.91	-21.0
-11			7	8	9	0	3	1	4.71	-25.5
-12			7	8	9	0	3	0	4.54	-27.2



Table 4. Computed pI's of some known proteins related to measured CPK pI's

	Protein Name	PIR Name	#ASP 3.9	#GLU 4.1	#HIS 6.0	#LYS 10.8	#ARG 12.5	Calc pI	R al CPK
0	Creatine phospho kinase (CPK), rabbit muscle	KIRBCM	28	27	17	34	18	6.84	0.0
1	Fatty acid-binding protein, rat hepatic	FZRTL	5	13	2	16	2	7.83	-3.0
2	b2-microglobulin, human	MGHUB2	7	8	4	8	5	6.09	-5.0
3	Carbamoyl-phosphate synthase, rat	SYRTCA	72	96	28	95	56	5.97	-5.5
4	Proalbumin ( serum albumin precursor), rat	ABRTS	32	57	15	53	27	5.98	-6.2
5	Serum albumin, rat	ABRTS	32	57	15	53	24	5.71	-9.0
6	Superoxid dismutase (Cu-Zn, SOD), rat	A26810	8	11	10	9	4	5.91	-9.2
7	Phospholipase C, phosphoinositide-specific (?), rat	A28807	34	42	9	49	21	5.92	-9.2
8	Albumin, human	ABHUS	36	61	16	60	24	5.70	-11.9
9	Apo A-I lipoprotein, rat	A24700	18	24	6	23	12	5.32	-13.7
10	proApo A-I lipoprotein, human	LPHUA1	16	30	6	21	17	5.35	-14.3
11	NADPH cytochrome P-450 reductase, rat	RDRT04	41	60	21	38	36	5.07	-15.6
12	Retinol binding protein, human	VAHU	18	10	2	10	14	5.04	-16.9
13	Actin beta, rat	ATRTC	23	26	9	19	18	5.06	-17.2
14	Actin gamma, rat	ATRTC	20	29	9	19	18	5.07	-16.8
15	Apo A-I lipoprotein, human	LPHUA1	16	30	5	21	16	5.10	-17.5
16	Apo A-IV lipoprotein, human	LPHUA4	20	49	8	28	24	4.88	-19.7
17	Tubulin alpha, rat	UBRTA	27	37	13	19	21	4.66	-19.8
18	F1ATPase beta, bovine	PWBOB	25	36	9	22	22	4.80	-21.0
19	Tubulin beta, pig	UBPGB	26	36	10	15	22	4.49	-22.5
20	Protein disulphide isomerase (PDI), rat hepatic	ISRTSS	43	51	11	51	9	4.07	-25.0
21	Cytochrome b5, rat	CBRT5	10	15	6	10	4	4.59	-26.0
22	Apo C-II lipoprotein, human	LPHUC2	4	7	0	6	1	4.44	-30.5
Amino acid pI assumed in calculation:			3.9	4.1	6.0	10.8	12.5		

# ELECTROPHORESIS

An International Journal

1191

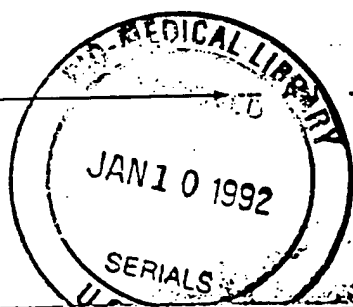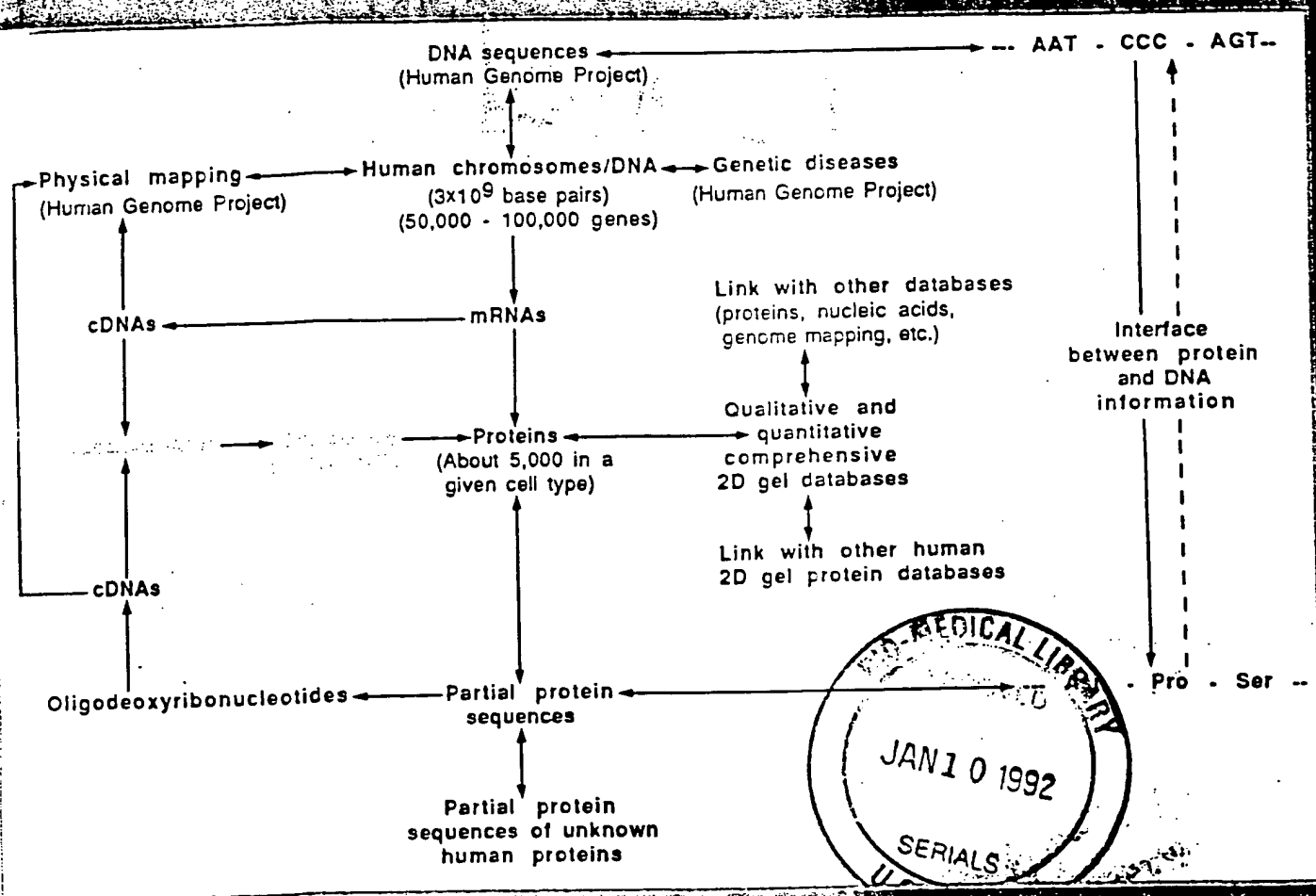

## TWO DIMENSIONAL GEL PROTEIN DATABASES

Editor: J. E. Celis



# ELECTROPHORESIS

An International Journal

Indexed in: BIOSIS  
Current Contents, MEDLAR  
ISSN 0173-0835  
ELCTDN 12 (11) 763-996 (1991)

## TWO-DIMENSIONAL GEL PROTEIN DATABASES

Editor: J. E. Celis

### Editorial

- J. E. Celis, H. Leffers,  
H. H. Rasmussen, P. Madsen,  
B. Honoré, B. Gesser, K. Dejgaard,  
E. Olsen, G. P. Ratz, J. B. Lauridsen,  
B. Basse, A. H. Andersen,  
E. Walbum, B. Brandstrup, A. Celis,  
M. Puype, J. Van Damme and  
J. Vandekerckhove
- 765 The master two-dimensional gel database of human AMA cell proteins: Towards linking protein and genome sequence and mapping information (Update 1991)
- J. E. Celis, P. Madsen,  
H. H. Rasmussen, H. Leffers,  
B. Honoré, B. Gesser, K. Dejgaard,  
E. Olsen, N. Magnusson, J. Kill,  
A. Celis, J. B. Lauridsen, B. Basse,  
G. P. Ratz, A. H. Andersen,  
E. Walbum, B. Brandstrup,  
P. S. Pedersen, N. J. Brandt,  
M. Puype, J. Van Damme and  
J. Vandekerckhove
- 802 A comprehensive two-dimensional gel protein database of noncultured unfractionated normal human epidermal keratinocytes: Towards an integrated approach to the study of cell proliferation, differentiation and skin diseases
- H. H. Rasmussen, J. Van Damme,  
M. Puype, B. Gesser, J. E. Celis and  
J. Vandekerckhove
- 873 Microsequencing of proteins recorded in human two-dimensional gel protein databases
- N. L. Anderson and N. G. Anderson
- 883 A two-dimensional gel database of human plasma proteins
- N. L. Anderson, R. Esquer-Blasco,  
J.-P. Hofmann and N. G. Anderson
- 907 A two-dimensional gel database of rat liver proteins useful in gene regulation and drug effects studies
- P. J. Wirth, L.-di Luo, Y. Fujimoto,  
H. C. Bisgaard and A. D. Olson
- 931 The rat liver epithelial (RLE) cell protein database
- R. A. VanBogelen and  
F. C. Neidhardt
- 955 The gene-protein database of *Escherichia coli*: Edition 4
- 995 Miscellaneous

For submission of papers, see Instructions to Authors (last page of this issue)